


Polar Basics -- Beta Blockers at pH 10.5



15031

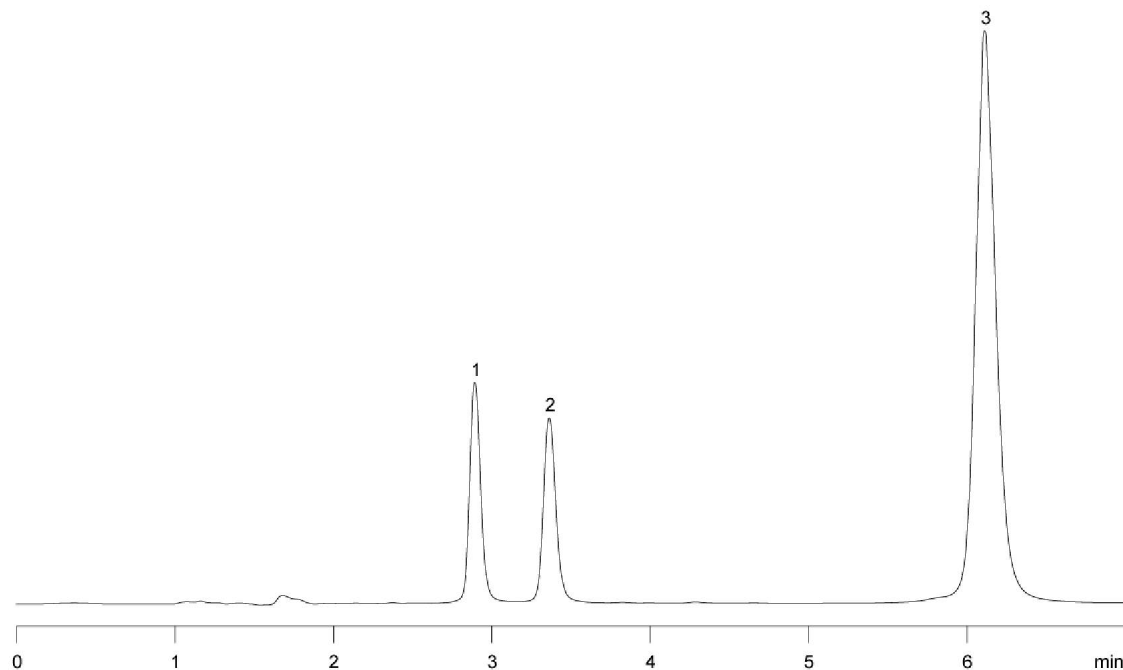

Column: Gemini 5 μ m C18 / ODS C18

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic

Eluent A: 10mM Ammonium Bicarbonate, pH10.5

Pct. A: 50

Eluent B: ACN

Pct. B: 50

Flow Rate: 1 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 230 nm

SAMPLE ANALYTES:

- 1 *Pindolol*
- 2 *Metoprolol*
- 3 *Propranolol*

SPE:

For more information call your Phenomenex Representative at

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 15031

Page 1 of 1