

Sorbitol-Mannitol-Maltitol on Luna Amino

Page 1 of 1

Column: Luna 5u NH2 100A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4378-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Pct. A: 85


Eluent B: Water

Pct. B: 15

Flow Rate: 3 mL/min


Col. Temp.: 40 °C

Detection: Refractive Index (RI) (ambient)

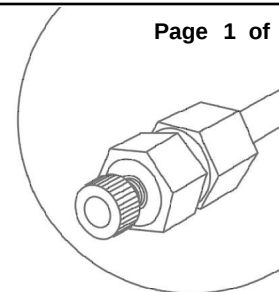

Products used in this Application:

LUNA®

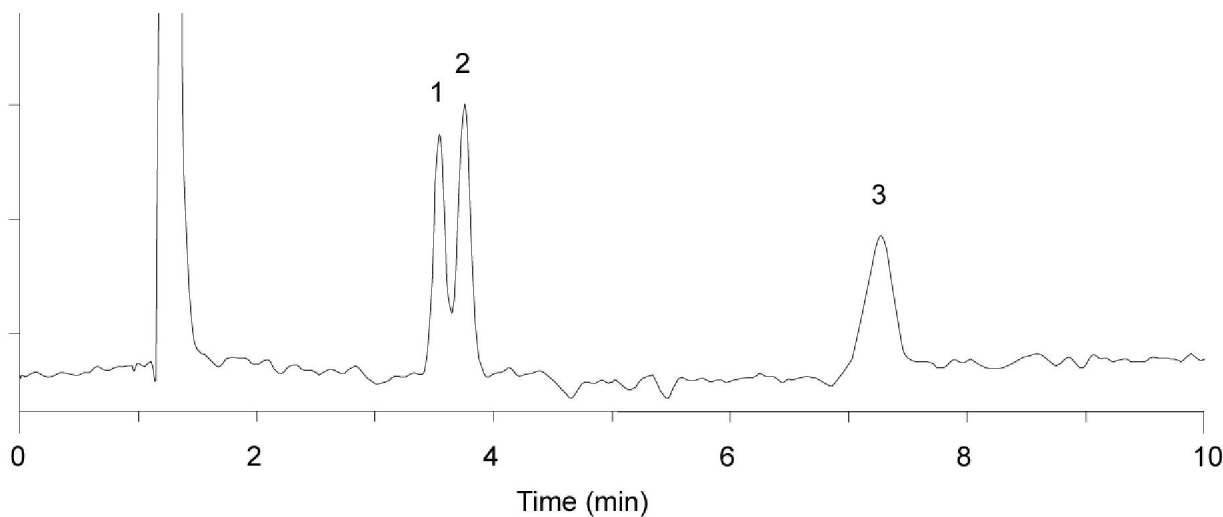

ANALYTES:

- 1 Sorbitol
- 2 Mannitol
- 3 Maltitol



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.