


Green Tea

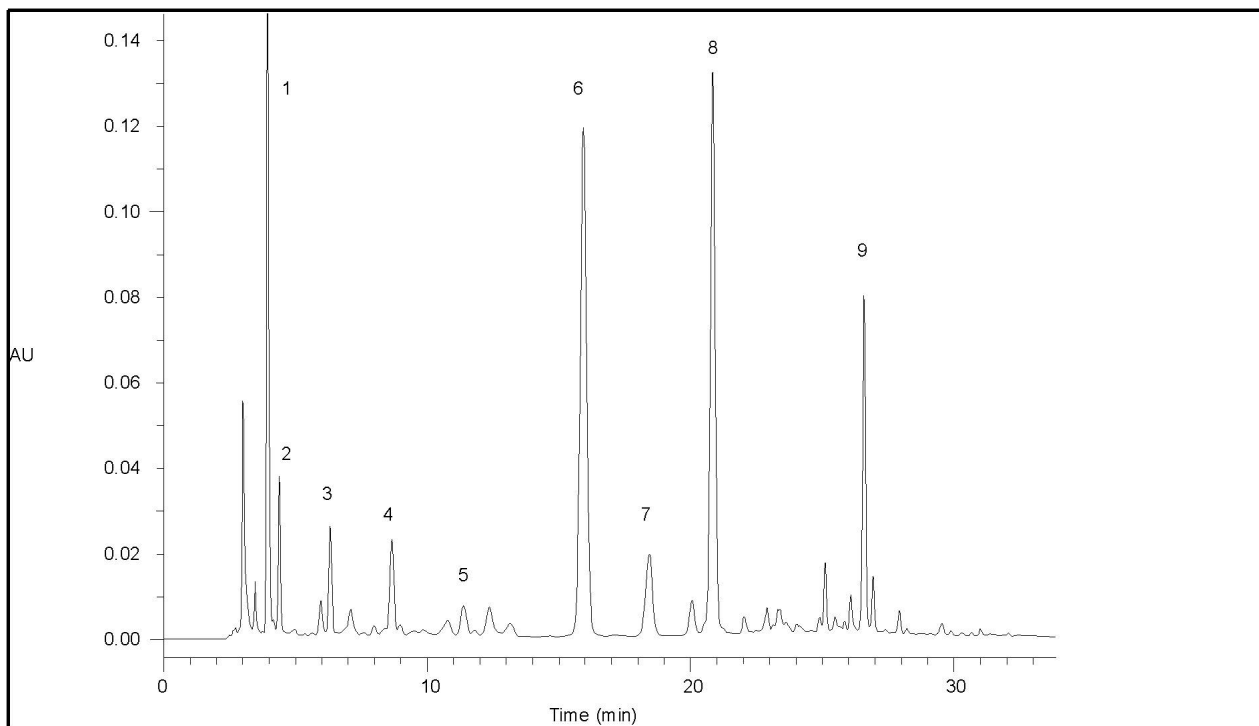

Column: Luna 5 μ m C6H5 / Phenyl Phenyl-Hexyl
Dimensions: 250 x 4.6 mm ID
Order No:
Elution Type: Gradient
Eluent A: 0.2% Acetic acid in Water/Acetonitrile (91:9)
Eluent B: Water/Acetonitrile (20:80)

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	100	0
2	10	100	0
3	25	68	32
4	35	68	32
5	40	100	0
6	60	100	0

Flow Rate: 1 mL/min
Temperature: 35 °C
Detection: UV-Vis Abs. - Variable Wave. @ 278 nm

SPE:

SAMPLE ANALYTE

1 Theogalline	5 Catechine
2 Gallic acid	6 Caffeine
3 Theobromine	7 Epicatechine
4 Epigallocatechine	8 Epigallocatechingallate


For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14436

Continued...

Green Tea

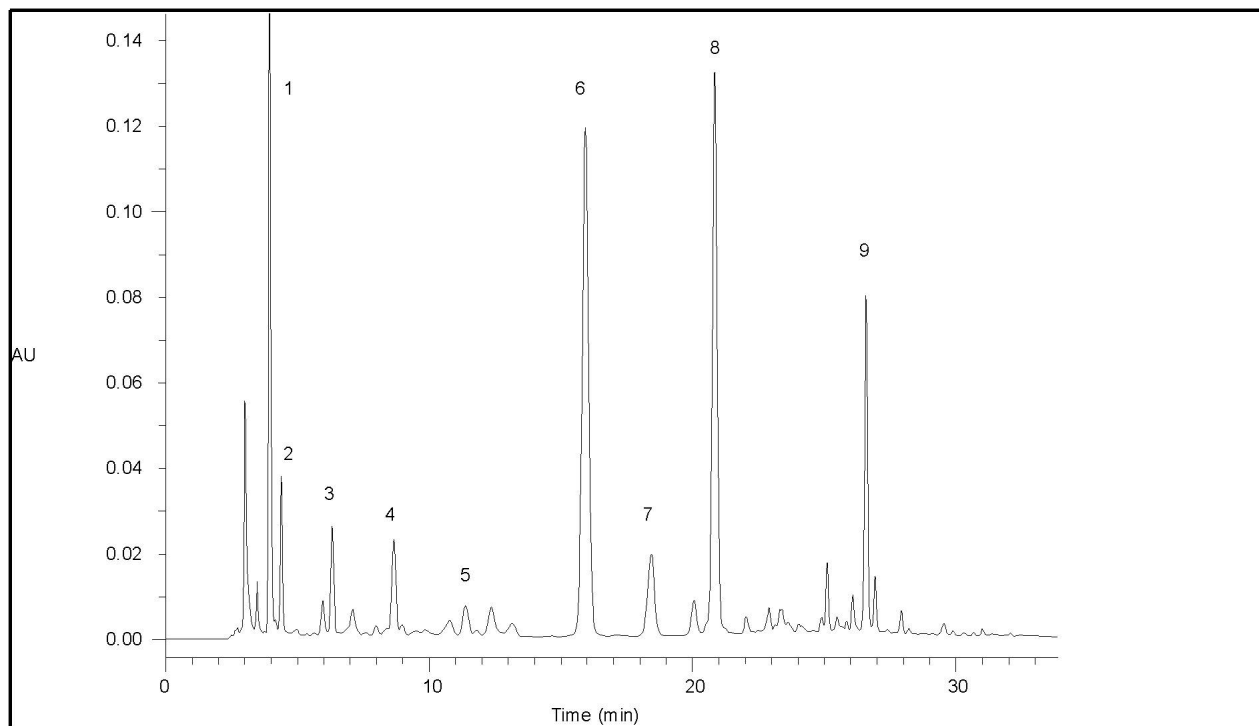

9 Epicatechingallate

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14436