


EPA Method 8316, Acrylamide, Acrylonitrile, & Acrolein

Page 1 of 1

Column: Synergi 4u Hydro-RP 80A
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4375-E0
Elution Type: Isocratic
Eluent A: Water
Eluent B: Acetonitrile
Flow Rate: 2 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs. - Variable Wave. @ 195 nm (ambient)


Pct. A: 98

Pct. B: 2

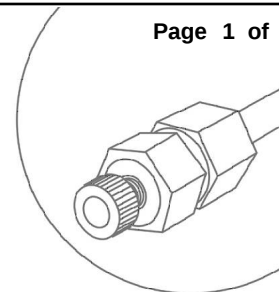

Products used in this Application:

synergi™

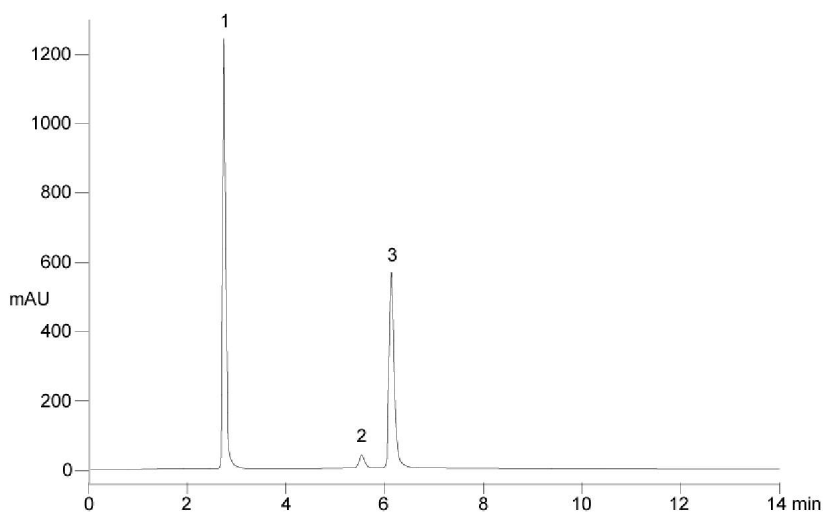

ANALYTES:

- 1 Acrylamide
- 2 Acrylonitrile
- 3 Acrolein



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.