


Gemini - High pH Prep Loading in DMSO



15464

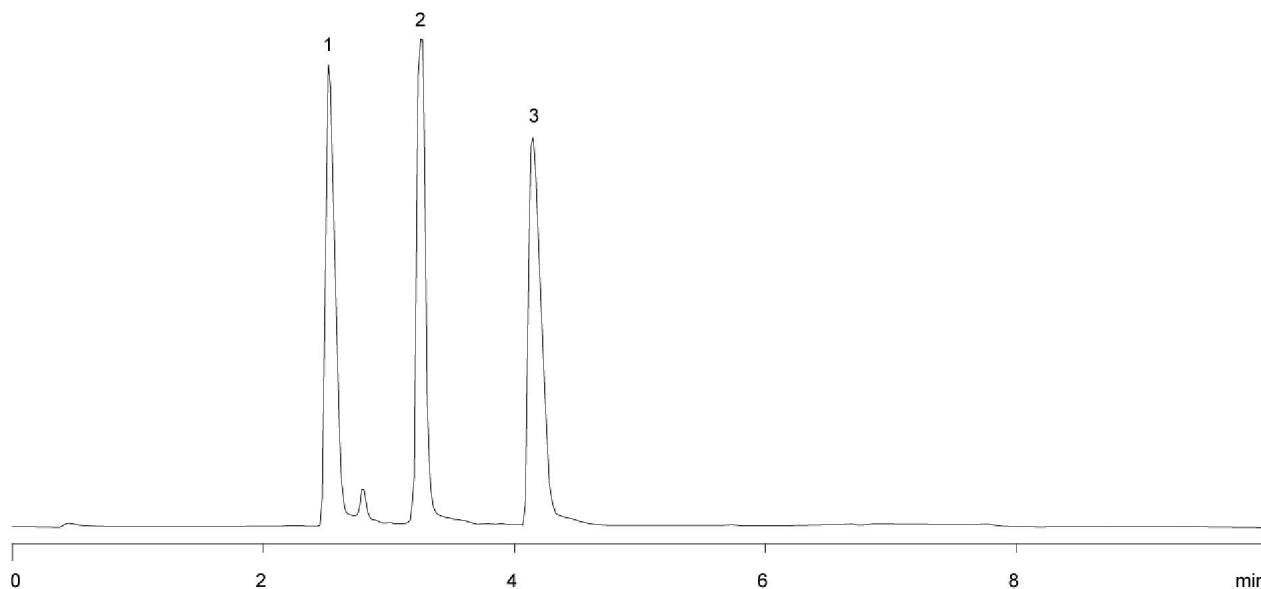

Column: Gemini 5 μ m C18 / ODS C18
Dimensions: 50 x 21.2 mm ID
Order No: 00B-4435-P0
Elution Type: Gradient
Eluent A: 90:10 50mM Ammonium Bicarbonate, pH 10/Acetonitrile
Eluent B: 10:90 50mM Ammonium Bicarbonate, pH 10/Acetonitrile

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	100	0
2	7	0	100
3	7.01	100	0
4	10	100	0

Flow Rate: 34 mL/min
Temperature: 22 °C
Detection: UV-Vis Abs. - Variable Wave. @ 254 nm
SPE:

SAMPLE ANALYTES:

- 1 Diphenhydramine
- 2 Oxybutynin chloride
- 3 Terfenadine

For more information call your Phenomenex Representative at

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 15464

Page 1 of 1