

SFC Application of Luna 10u CN

Column: Luna 10u CN 100A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4300-E0

Elution Type: Isocratic

Eluent A: Supercritical Carbon Dioxide

Eluent B: Methanol

Flow Rate: 4 mL/min


Col. Temp.: 40 °C

Detection: UV-Vis Abs. - Variable Wave. @ 265 nm (40 °C)

Pct. A: 85


Pct. B: 15

Page 1 of 1

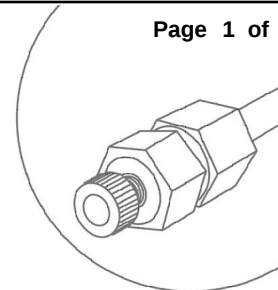

Products used in this Application:

LUNA™

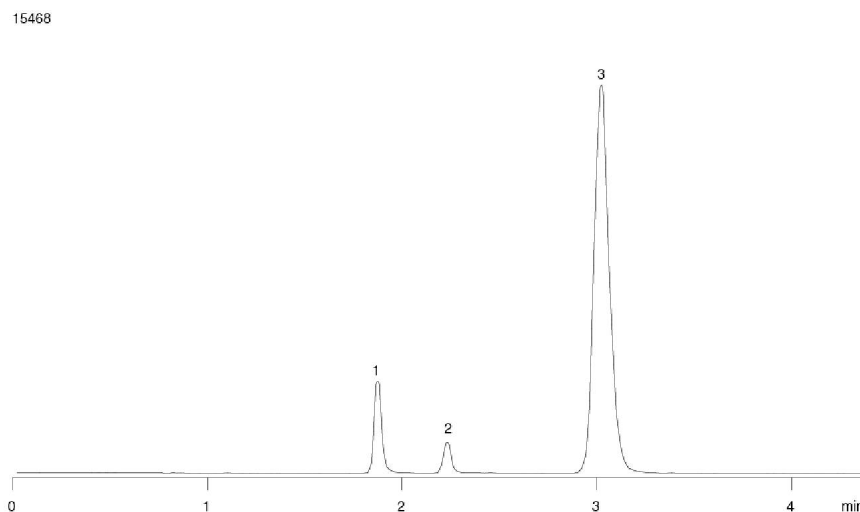

ANALYTES:

- 1 Acetaminophen
- 2 Benzoic Acid
- 3 Procaine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.