

Beta Blockers on Onyx C18 - Effect of Flow Rate

Column: Onyx Monolithic C18

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7643

Elution Type: Isocratic

Eluent A: 0.1% TFA in Water

Eluent B: Acetonitrile

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 220 nm (25 °C)

Pct. A: 80


Pct. B: 20

Page 1 of 1

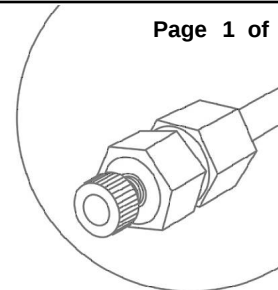

Products used in this Application:

ONYX

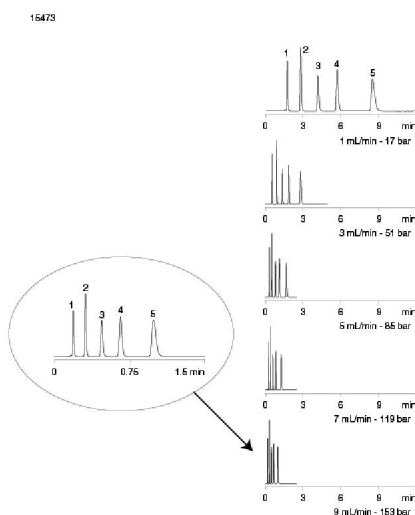

ANALYTES:

- 1 Atenolol
- 2 Pindolol
- 3 Metoprolol
- 4 Celiprolol
- 5 Bisoprolol

AUTHORIZED DISTRIBUTOR



SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr