

Beta Blockers on Onyx C18

Column: Onyx Monolithic C18

Dimensions: 50 x 4.6 mm ID

Order No: CH0-7644

Elution Type: Isocratic

Eluent A: 0.1% TFA in Water

Pct. A: 80

Eluent B: Acetonitrile


Pct. B: 20

Flow Rate: 5 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 220 nm (ambient)

Page 1 of 1

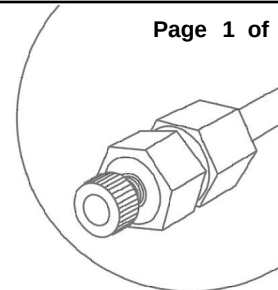

Products used in this Application:

ONYX

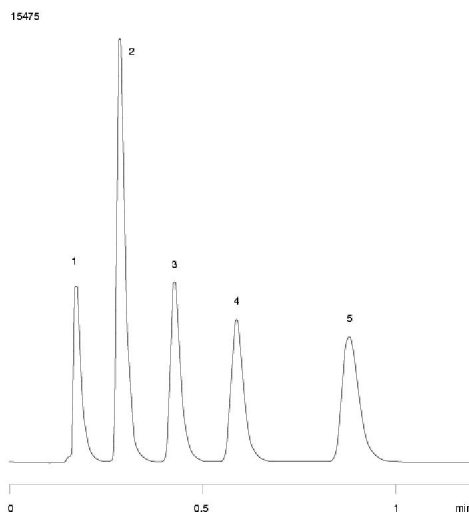

ANALYTES:

- 1 Atenolol
- 2 Pindolol
- 3 Metoprolol
- 4 Celiprolol
- 5 Bisoprolol

AUTHORIZED DISTRIBUTOR



SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr