

Phthalates on Onyx C18

Column: Onyx Monolithic C18

Dimensions: 50 x 4.6 mm ID

Order No: CH0-7644

Elution Type: Gradient

Eluent A: Water


Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40
	2	0.1	0	100
	3	2	0	100


Flow Rate: 4 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

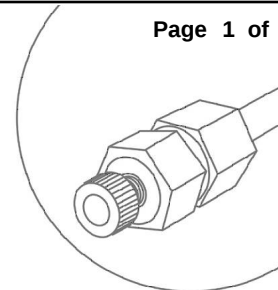

Products used in this Application:

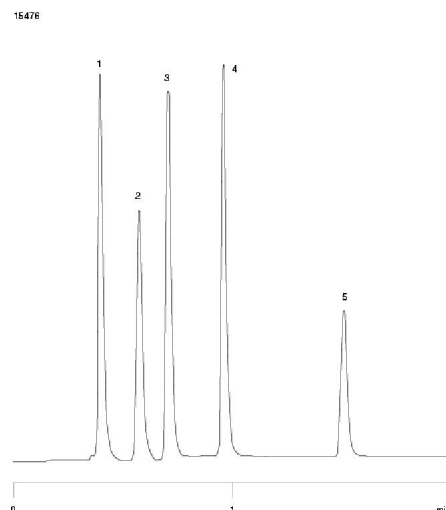ONYX

ANALYTES:

- 1 Dimethyl Phthalate
- 2 Dimethyl Terephthalate
- 3 Diethyl Phthalate
- 4 Benzyl Butyl Phthalate
- 5 Bis(2-ethylhexyl) Phthalate

AUTHORIZED DISTRIBUTOR



SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr