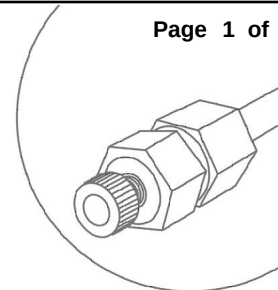
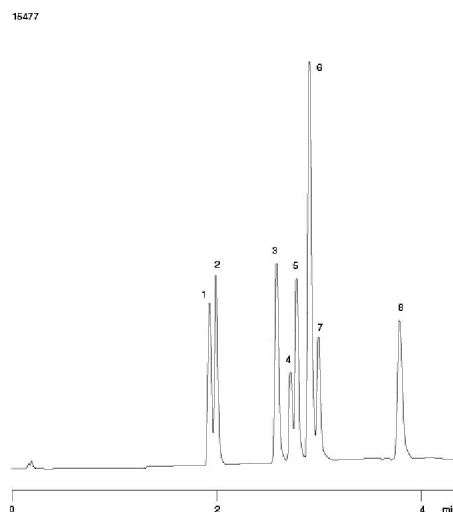


Steroids on Onyx C18**Column:** Onyx Monolithic C18**Dimensions:** 50 x 4.6 mm ID**Order No:** CH0-7644**Elution Type:** Gradient**Eluent A:** Water**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	3	50	50
	3	4.5	50	50

Flow Rate: 4 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 220 nm (ambient)

Products used in this Application:

ONYX**ANALYTES:**

- 1 Prednisolone
- 2 Cortisone
- 3 Nortestosterone
- 4 Estradiol
- 5 Testosterone
- 6 Corticosterone
- 7 Estrone
- 8 Progesterone

AUTHORIZED DISTRIBUTOR

SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr