

Cationic Surfactants on Onyx C18

Column: Onyx Monolithic C18

Dimensions: 50 x 4.6 mm ID

Order No: CH0-7644

Elution Type: Gradient

Eluent A: 0.1% TFA in Water

Eluent B: Methanol


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	40	60
	2	2.5	0	100

Flow Rate: 4 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

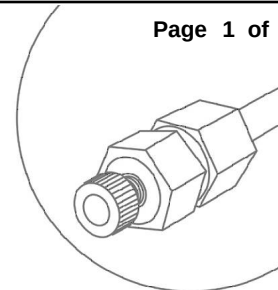

Products used in this Application:

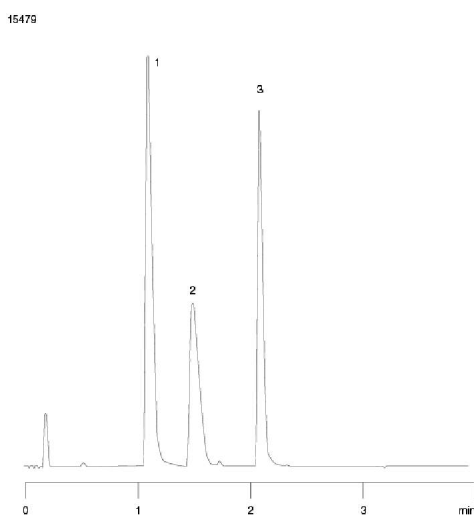ONYX[™]

ANALYTES:

- 1 N-Dodecylpyridinium Chloride
- 2 Hyamine
- 3 N-Cetylpyridinium Chloride

AUTHORIZED DISTRIBUTOR



SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr