

Peptides (Angiotensin) on Onyx C8

Page 1 of 1

Column: Onyx Monolithic C8

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7647

Elution Type: Gradient

Eluent A: 0.1% TFA in Water

Eluent B: 0.1% TFA in 90/10 Acetonitrile:Water

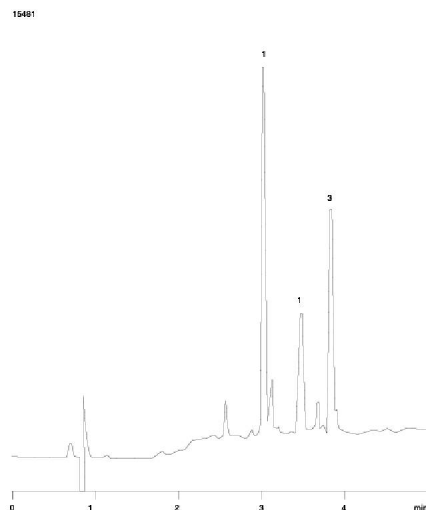
Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	55	45
	2	1	10	90
	3	3	10	90

Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 214 nm (ambient)

Products used in this Application:

ONYX[™]

ANALYTES:

- 1 (Sar1, Ala8) Angiotensin II
- 2 (Sar1, Ile8) Angiotensin II
- 3 Angiotensin



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.