

Chelators on Onyx C8

Column: Onyx Monolithic C8

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7647

Elution Type: Isocratic

Eluent A: 20mM Phosphate buffer, pH 2.0

Eluent B: Methanol

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Pct. A: 30


Pct. B: 70

Page 1 of 1

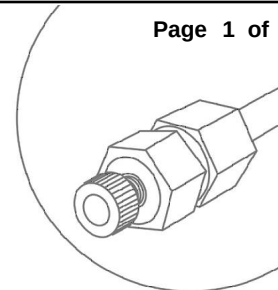

Products used in this Application:

ONYX

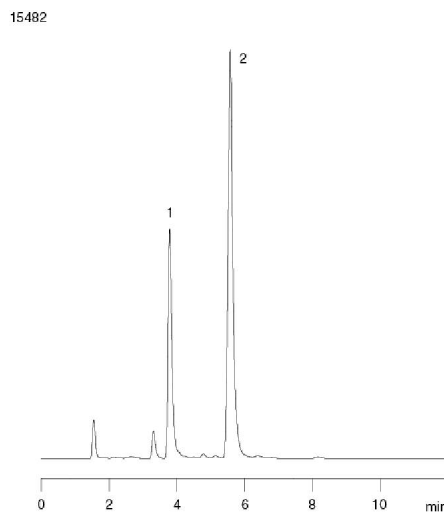

ANALYTES:

- 1 Purpurin
- 2 Quinizarin

AUTHORIZED DISTRIBUTOR



SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr