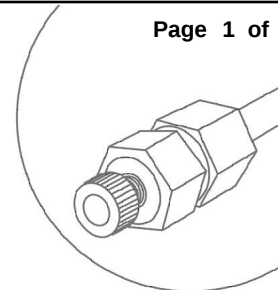
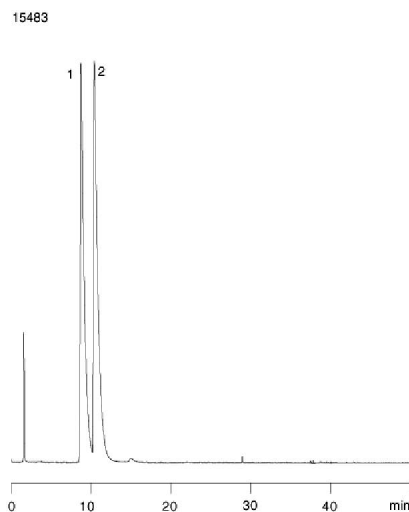


Triptylines on Onyx C8

Column: Onyx Monolithic C8**Dimensions:** 100 x 4.6 mm ID**Order No:** CH0-7647**Elution Type:** Isocratic**Eluent A:** 20mM Phosphate buffer, pH 7.6**Pct. A:** 30**Eluent B:** Methanol**Pct. B:** 70**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

**Products used in this Application:****ONYX****ANALYTES:**

- 1 Imipramine
- 2 Amitriptyline

**www.phenomenex.com**

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.