

Anisoles on Onyx Silica - Effect of Flow Rate

Column: Onyx Monolithic Si

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7648

Elution Type: Isocratic

Eluent A: n-Heptane

Eluent B: Dioxane

Flow Rate: 2 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Pct. A: 95

Pct. B: 5





Page 1 of 1

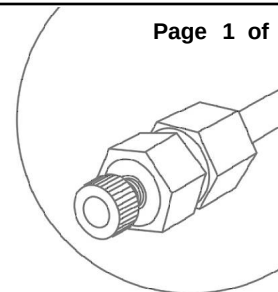

Products used in this Application:

ONYX

15484

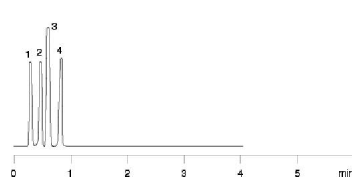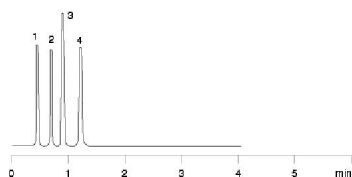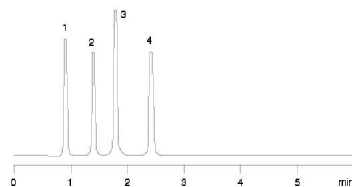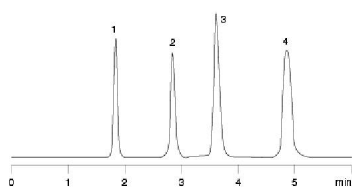

ANALYTES:

- 1 Anisole
- 2 3-Nitroanisole
- 3 4-Nitroanisole
- 4 2-Nitroanisole

AUTHORIZED DISTRIBUTOR



SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr