

Peptides from Serum on Strata™ X

Page 1 of 1

Column: Jupiter 4u Proteo 90A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4396-E0

Elution Type: Gradient

Eluent A: 0.1% TFA

Eluent B: ACN:0.1% TFA in water (90:10)

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	85	15
	2	22	40	60

Flow Rate: 1 mL/min

Col. Temp.: 30 °C

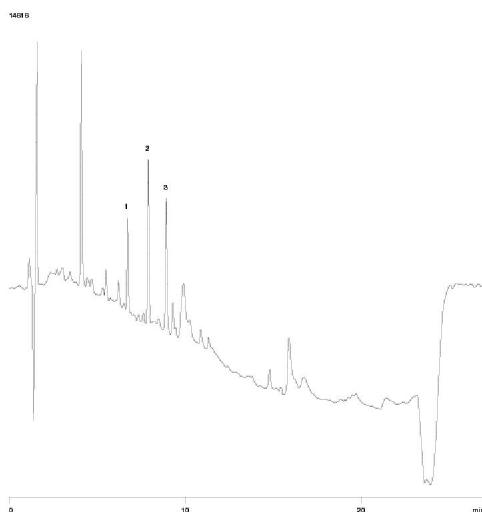
Detection: UV-Vis Abs. - Diode Array @ 214 nm (ambient)

Solid Phase Extraction: strata-X 33µm Polymeric Reversed Phase, 30 mg / 1 mL, Tube, 100/Pk

Order No.: 8B-S100-TAK

(SPE): See Sample Preparation report for full details of the extraction procedure

Products used in this Application:

JupiterstrataX

ANALYTES:

- 1 Oxytocin
- 2 Angiotensin II
- 3 Angiotensin I



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au