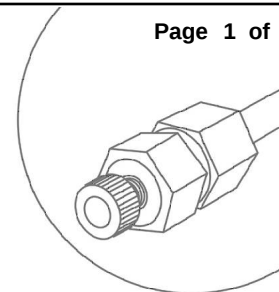


Hydrophobic Bases at Low pH on Onyx C18


Page 1 of 1

Column: Onyx Monolithic C18**Dimensions:** 100 x 4.6 mm ID**Order No:** CH0-7643**Elution Type:** Isocratic**Eluent A:** 20mM Phosphate buffer, pH 2.5**Pct. A:** 25**Eluent B:** Acetonitrile**Pct. B:** 75**Flow Rate:** 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Products used in this Application:

ONYX

15486

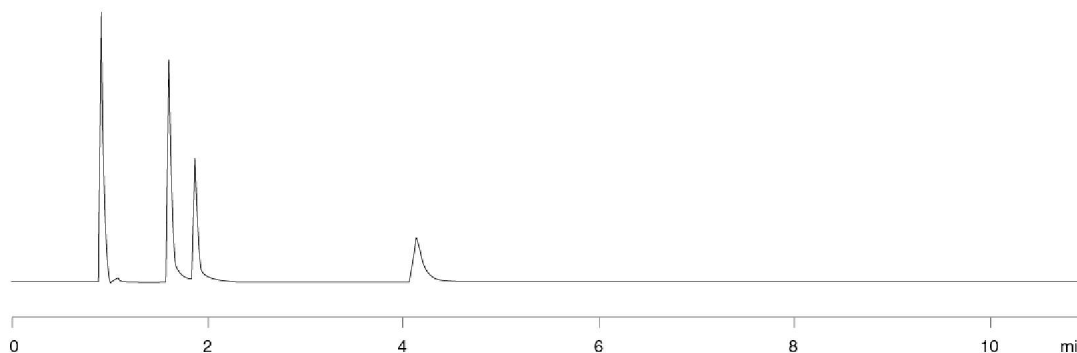**ANALYTES:**

- 1 Maleic Acid
- 2 Triprolidine
- 3 Chlorpheniramine
- 4 Diphenhydramine

AUTHORIZED DISTRIBUTOR

SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr