


Beta Blockers on Onyx C18

Column: Onyx Monolithic C18**Dimensions:** 100 x 4.6 mm ID**Order No:** CH0-7643**Elution Type:** Isocratic**Eluent A:** 20mM Sodium Phosphate, pH 3.0**Pct. A:** 70**Eluent B:** Methanol**Pct. B:** 30**Flow Rate:** 2 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)


Page 1 of 1

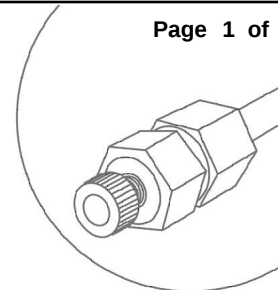

Products used in this Application:

ONYX

15487

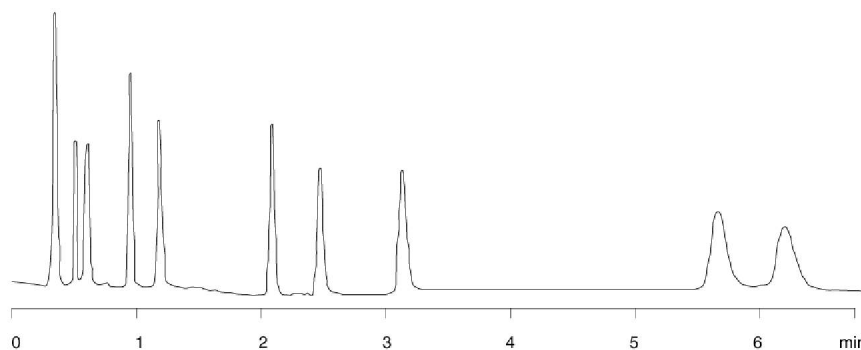

ANALYTES:

- 1 Fumaric acid
- 2 Pindolol
- 3 Nadolol
- 4 Pafenolol
- 5 Metoprolol
- 6 Celiprolol
- 7 Carazolol
- 8 Bisoprolol
- 9 Alprenolol
- 10 Propranolol

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.