

Phenols on Onyx C8

Page 1 of 1

Column: Onyx Monolithic C8

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7647

Elution Type: Gradient

Eluent A: 0.1% Phosphoric acid

Eluent B: Acetonitrile

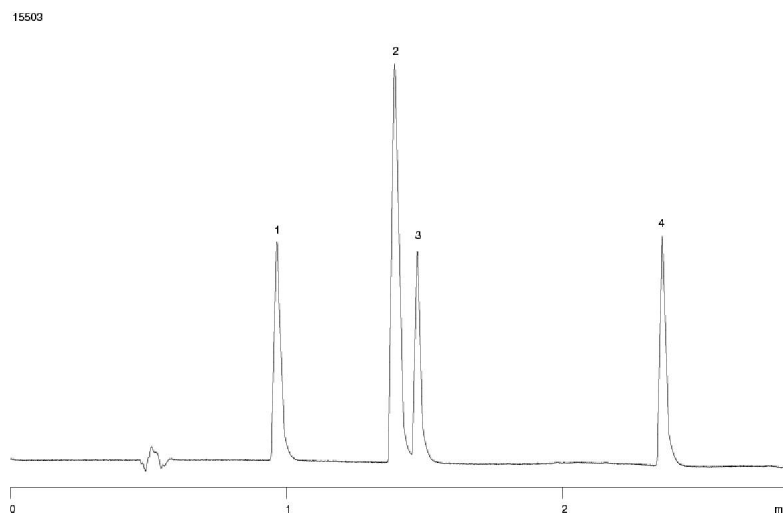
Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30
	2	2	20	80

Flow Rate: 3 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Products used in this Application:

ONYX[™]

ANALYTES:

- 1 Phenol
- 2 2,4-Dinitrophenol
- 3 3,4-Dinitrophenol
- 4 Pentachlorophenol

AUTHORIZED DISTRIBUTOR



SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr