

Vitamins on Onyx C8

Column: Onyx Monolithic C8

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7647

Elution Type: Gradient

Eluent A: 20mM Potassium phosphate, pH 3.0


Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	6	75	25

Flow Rate: 2 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

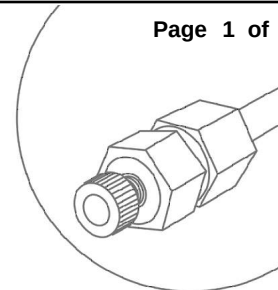

Products used in this Application:

ONYX

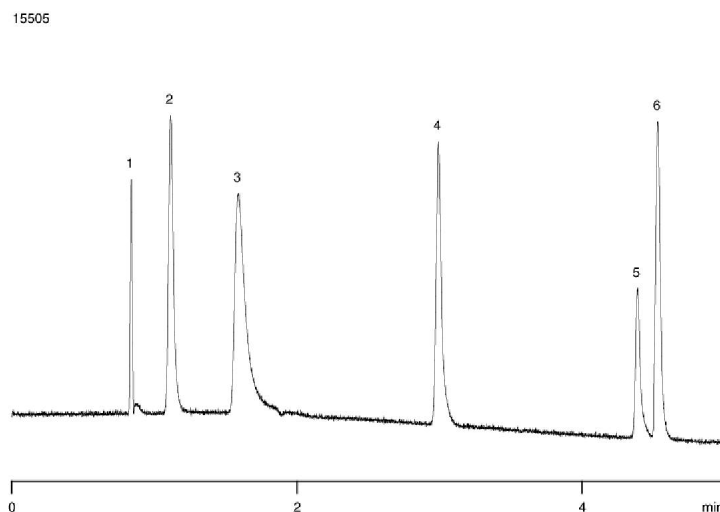

ANALYTES:

- 1 Thiamine
- 2 Uracil
- 3 Niacinamide
- 4 p-Aminobenzoic acid
- 5 Riboflavin (Vitamin B2)
- 6 Vitamin B12

