


Energy Drink Components on Onyx C18

Column: Onyx Monolithic C18**Dimensions:** 50 x 4.6 mm ID**Order No:** CH0-7644**Elution Type:** Isocratic**Eluent A:** 10mM Potassium phosphate, pH 5.0**Eluent B:** Acetonitrile**Flow Rate:** 4 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 227 nm (ambient)**Pct. A:** 97**Pct. B:** 3


Page 1 of 1

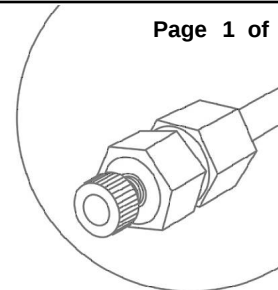

Products used in this Application:

ONYX

15506

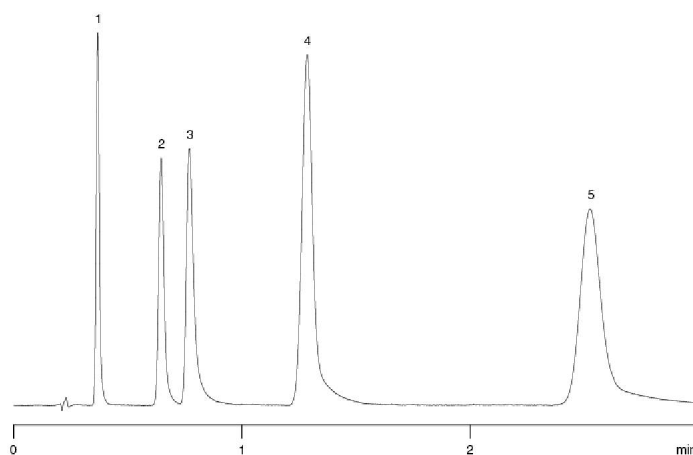**ANALYTES:**

- 1 p-Hydroxybenzoic acid
- 2 Saccharin
- 3 Salicylic acid
- 4 Sorbic acid
- 5 p-Toluic acid

AUTHORIZED DISTRIBUTOR

SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr