

Geometric Isomers on Onyx C18 @ 3mL/min

Column: Onyx Monolithic C18

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7643

Elution Type: Isocratic

Eluent A: 0.1%TFA in Water

Eluent B: ACN

Flow Rate: 3 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 220 nm (ambient)

Pct. A: 95

Pct. B: 5


Page 1 of 1

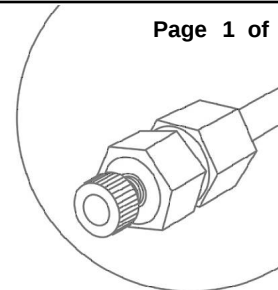

Products used in this Application:

ONYX

15507

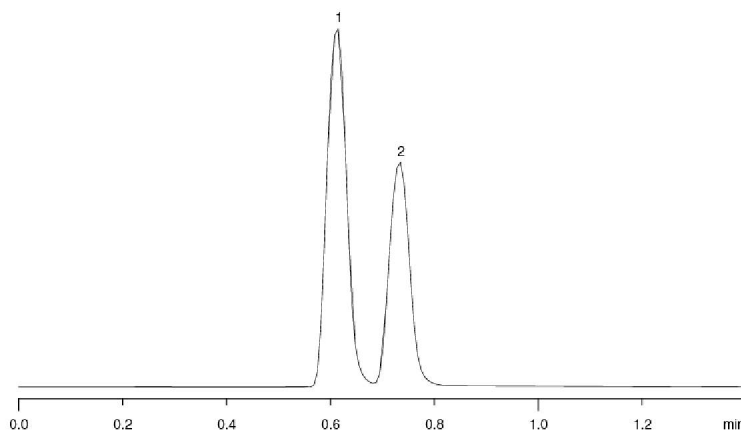

ANALYTES:

- 1 Maleic Acid
- 2 Fumaric Acid

AUTHORIZED DISTRIBUTOR



SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr