


30mer Oligo on Jupiter C18 300A

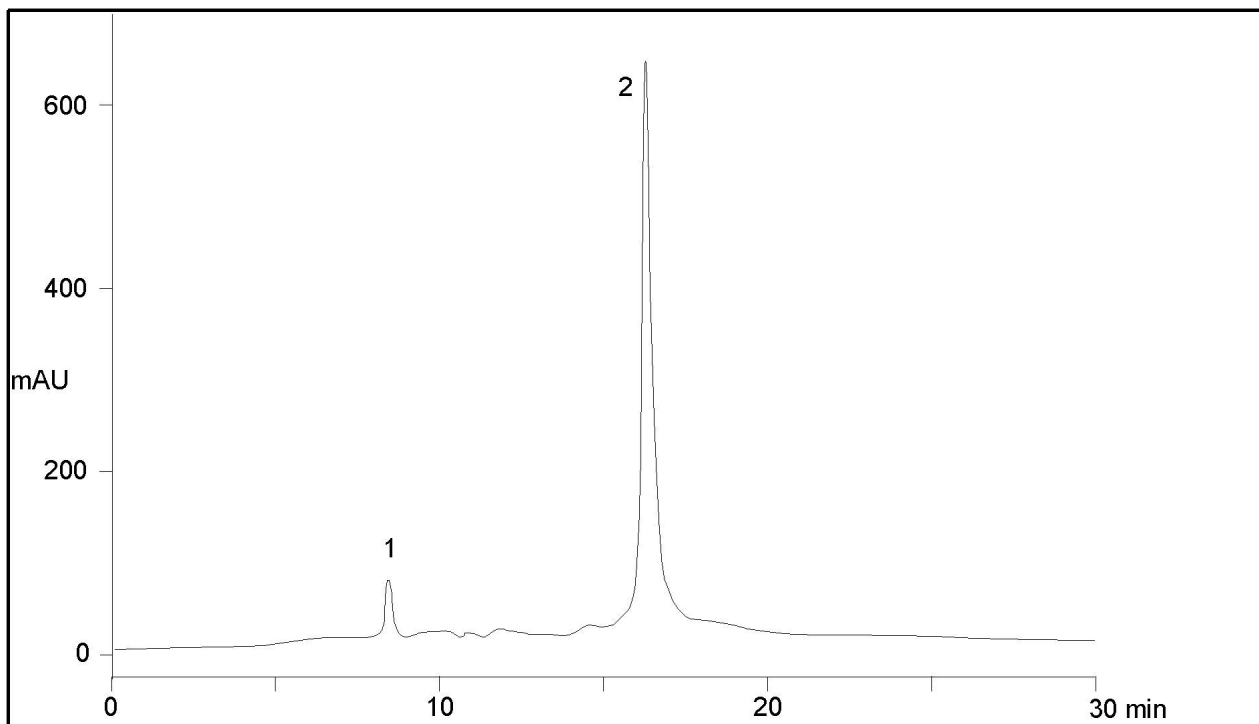

Column: Jupiter 5 μ m C18 / ODS
Dimensions: 150 x 4.6 mm ID
Order No:
Elution Type: Gradient
Eluent A: 100mM TEAA in Water, pH 7
Eluent B: 100mM TEAA in Acetonitrile, pH 7

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	1	96	4
2	30	85	15

Flow Rate: 1 mL/min
Temperature: 50 °C
Detection: UV-Vis Abs. - Diode Array @ nm
SPE:

SAMPLE ANALYTE

- 1 *Failed sequence*
- 2 *30mer Oligonucleotide*

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14624