

Pesticides on Onyx C18

Column: Onyx Monolithic C18

Dimensions: 300 x 4.6 mm ID

Order No:

Elution Type: Gradient

Eluent A: 10mM Ammonium acetate, pH 7.0

Eluent B: Acetonitrile


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	40	10	90

Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 214 nm (ambient)


Page 1 of 1

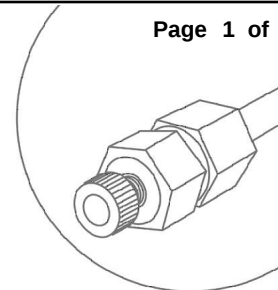

Products used in this Application:

ONYX[™]

Analyst Note: Catalog No. CH0-7643 (x 3), i.e., 3 columns in series

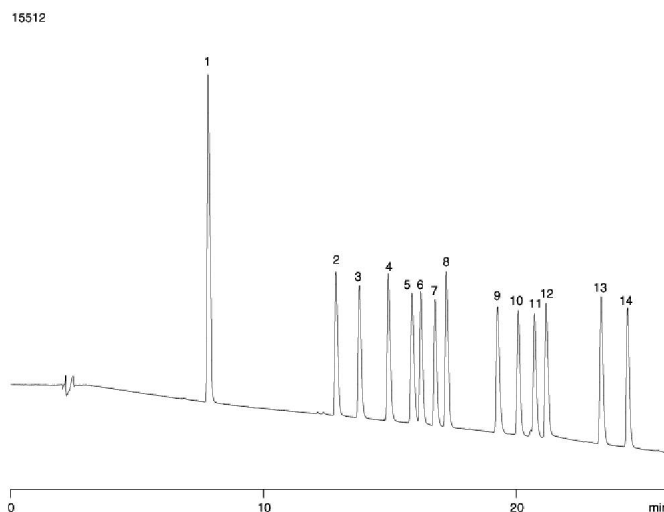

ANALYTES:

- 1 Fenuron
- 2 Monuron
- 3 Aminocarb
- 4 Carbofuran
- 5 Atrazine
- 6 Fluometuron
- 7 Forchlorfenuron
- 8 Diuron
- 9 Propazine
- 10 Siduron
- 11 Methiocarb
- 12 Linuron
- 13 Chlorpropham
- 14 Barban

©Copyright 2000-2005 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au