

Anisoles on Onyx Si

Column: Onyx Monolithic Si

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7648

Elution Type: Isocratic

Eluent A: Hexane

Pct. A: 95

Eluent B: Ethanol


Pct. B: 5

Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)


Page 1 of 1

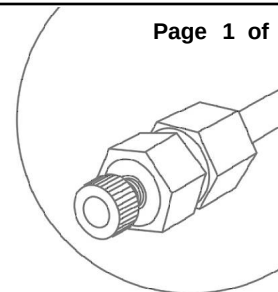

Products used in this Application:



15513

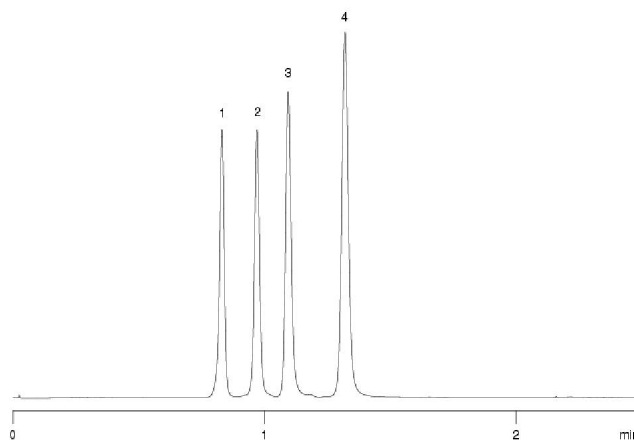

ANALYTES:

- 1 Anisole
- 2 3-Nitroanisole
- 3 4-Nitroanisole
- 4 2-Nitroanisole

AUTHORIZED DISTRIBUTOR



SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr