

Phenones on Onyx Si

Column: Onyx Monolithic Si

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7648

Elution Type: Isocratic

Eluent A: Hexane

Eluent B: Ethanol

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Pct. A: 95

Pct. B: 5


Page 1 of 1

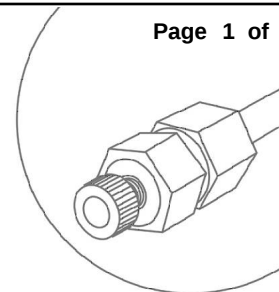

Products used in this Application:

ONYX

15514

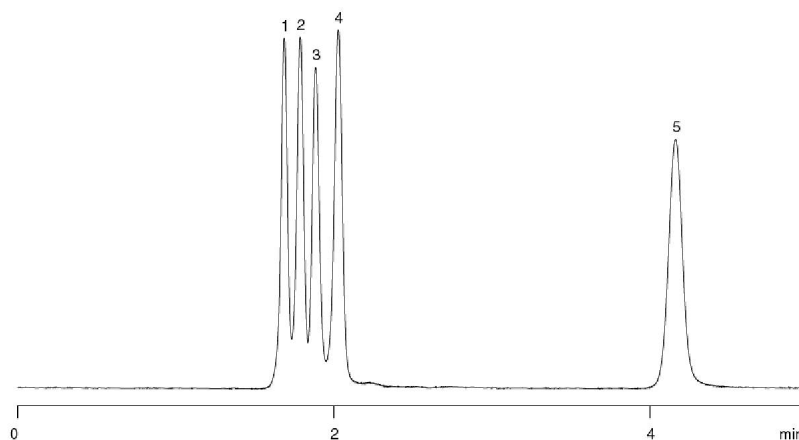

ANALYTES:

- 1 Octanophenone
- 2 Butyrophenone
- 3 Menadione
- 4 Acetophenone
- 5 Benzanilide

AUTHORIZED DISTRIBUTOR



SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr