

Anion Exchange 12-18mer Oligo on Luna Amino

Column: Luna 5u NH2 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4378-E0

Elution Type: Gradient

Eluent A: 20mM Sodium Phosphate in 90/10 Water/Acetonitrile

Eluent B: 20mM Sodium Phosphate in 90/10 Water/Acetonitrile + 1000mM NaCl


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	2	75	25
	2	40	5	95

Flow Rate: 1 mL/min

Col. Temp.: 40 °C


Detection: UV-Vis Abs. - Diode Array (ambient)

Page 1 of 1

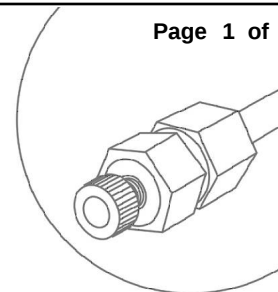

Products used in this Application:

LUNA™

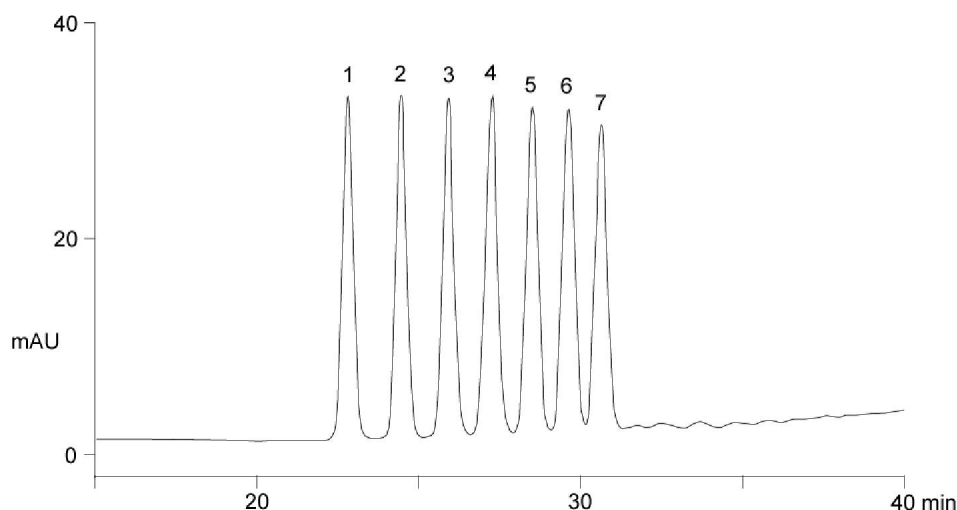

ANALYTES:

- 1 12mer Oligonucleotide
- 2 13mer Oligonucleotide
- 3 14mer Oligonucleotide
- 4 15mer Oligonucleotide
- 5 16mer Oligonucleotide
- 6 17mer Oligonucleotide
- 7 18mer Oligonucleotide

