

Separation of Aflatoxins using Gemini

Column: Gemini 5u C18 110A

Dimensions: 150 x 3 mm ID

Order No: 00F-4435-Y0

Elution Type: Isocratic

Eluent A: 10 mM Ammonium Acetate, pH 4.5

Eluent B: Acetonitrile

Flow Rate: 0.5 mL/min


Col. Temp.: ambient

Detection: Mass Spectrometer (MS) (ambient)

Pct. A: 70

Pct. B: 30


Page 1 of 1

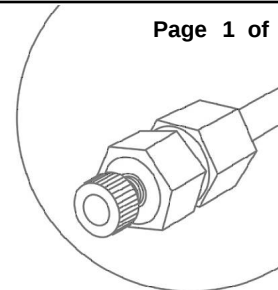

Products used in this Application:



15524

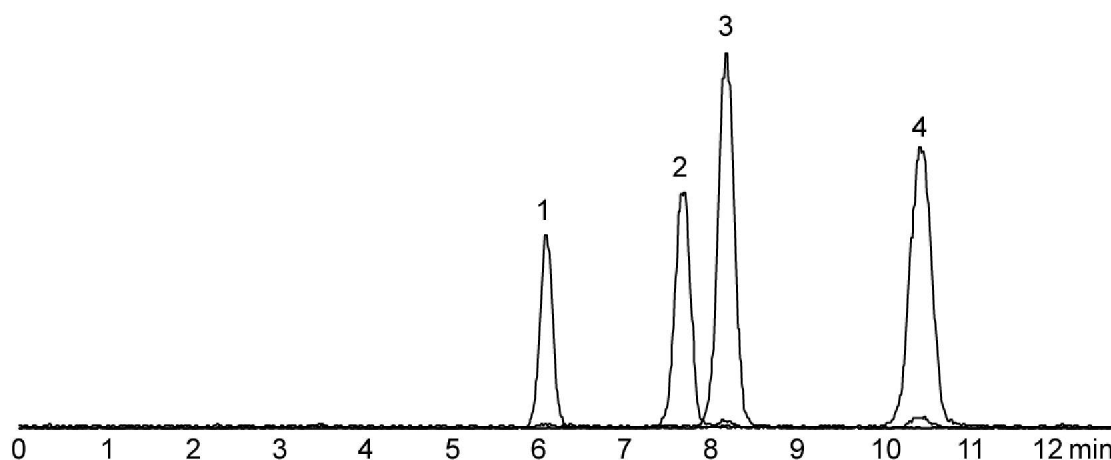

ANALYTES:

- 1 Aflatoxin G2
- 2 Aflatoxin B2
- 3 Aflatoxin G1
- 4 Aflatoxin B1

