


APPLICATION NOTE - HPLC
Polypeptides on PolymerX

polymerX™

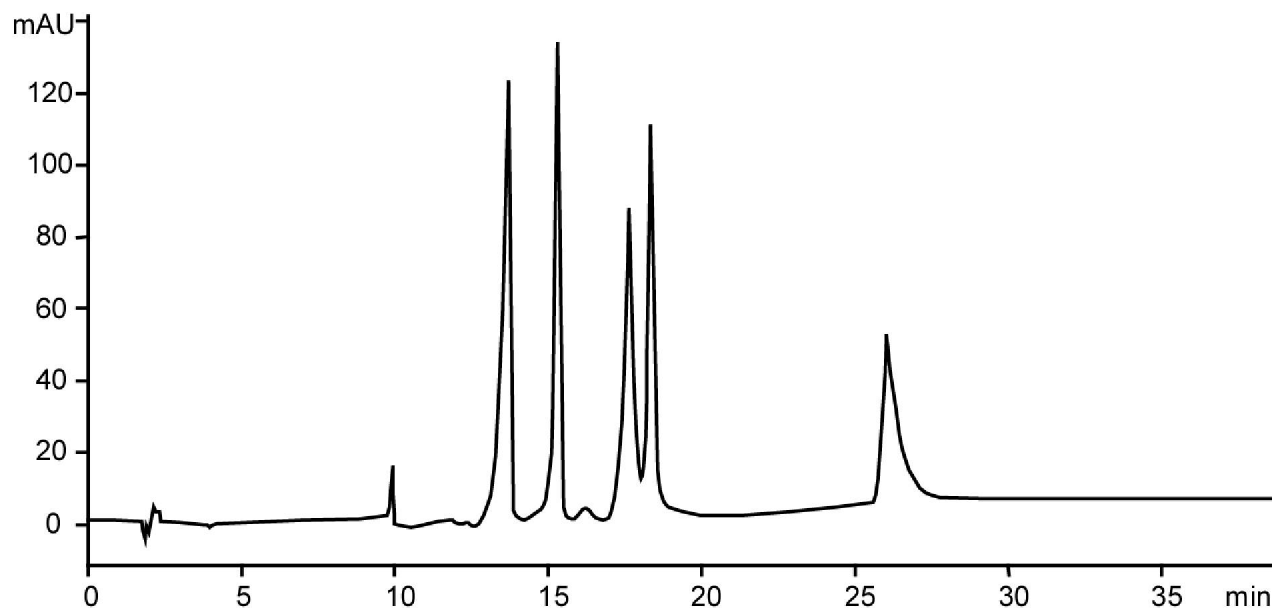

Column: PolymerX 5 μ m Unmodified copolymers of styrene divinylbenzene
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4326-E0
Elution Type: Gradient
Eluent A: 10mM NH₄OH
Eluent B: 20%: 10mM NH₄OH :80% ACN

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	90	10
2	25	55	45

Flow Rate: 0.8 mL/min
Temperature: 22 °C
Detection: UV-Vis Abs. - Variable Wave. @ nm
SPE:

SAMPLE ANALYTES:

- 1 *Bradykinin Frag 2-7*
- 2 *Bradykinin Frag 1-7*
- 3 *Bradyknin*
- 4 *Angiotensin I*
- 5 *Angiotensin II*

For more information call your Phenomenex Representative at (310) 212-0555

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 14762

Page 1 of 1