



HMW BioSep4000

Column: BioSep-SEC-S 4000**Dimensions:** 300 x 7.8 mm ID**Order No:** 00H-2147-K0**Elution Type:** Isocratic**Eluent A:** 20mM Na₂HPO₄**Pct. A:** 100**Flow Rate:** 0.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 280 nm (ambient)

Page 1 of 1

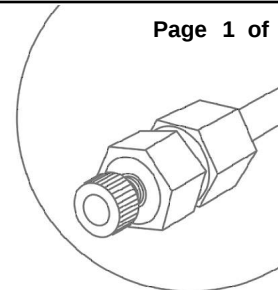

Products used in this Application:

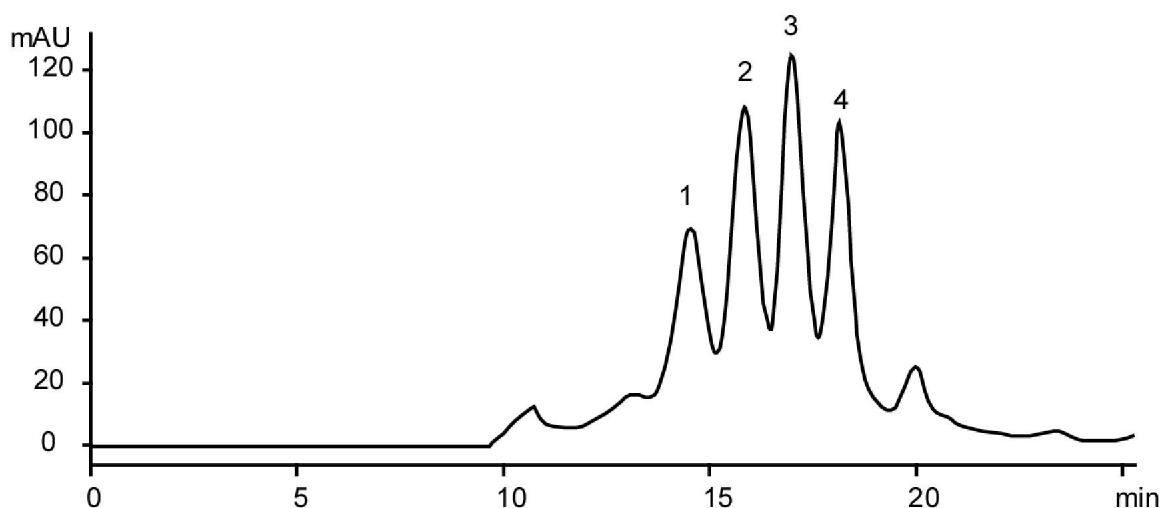Bio.Sep[™]**ANALYTES:**

- 1 Thyroglobulin
- 2 Apoferritin
- 3 Amylase
- 4 BSA

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.