

Trichothecenes analysis using Synergi Fusion-RP

Page 1 of 1

Column: Synergi 4u Fusion-RP 80A
Dimensions: 150 x 3 mm ID
Order No: 00F-4424-Y0
Elution Type: Gradient
Eluent A: 20mM Ammonium Acetate, pH 4.5
Eluent B: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0.5	90	10
	2	15	20	80

Flow Rate: 0.5 mL/min


Col. Temp.: ambient

Detection: Mass Spectrometer (MS) (ambient)

Products used in this Application:

synergi[™]

15541

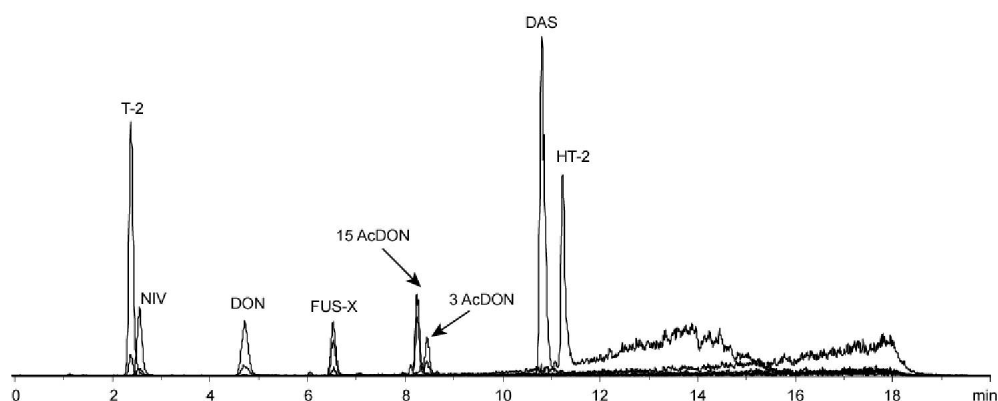

ANALYTES:

- 1 T-2
- 2 Nivalenol
- 3 Deoxynivalenol
- 4 Fusarenon-X
- 5 15-Acety-Deoxynivalenol
- 6 3-Acetyl-Deoxynivalenol
- 7 Diacetoxyscirpenol
- 8 HT-2



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au