

siRNA on Jupiter Proteo

Column: Jupiter 4u Proteo 90A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4396-E0

Elution Type: Gradient

Eluent A: 100mM TEAA pH 7


Eluent B: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0.5	96	4
	2	25	88	12


Flow Rate: 1 mL/min

Col. Temp.: 30 °C

Detection: UV-Vis Abs. - Variable Wave. (ambient)

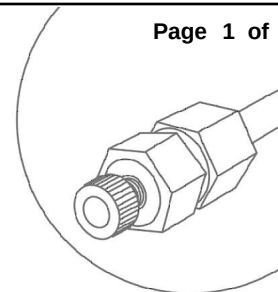

Products used in this Application:

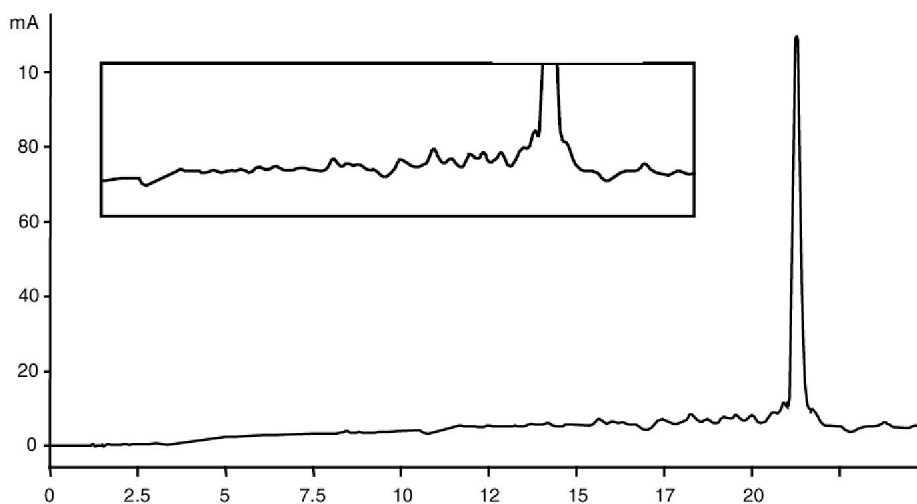

ANALYTES:

1 siRNA 19mer



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.