


Beta Amylase Digest**Column:** Onyx Monolithic C18**Dimensions:** 150 x 0.1 mm ID**Order No:** CH0-7646**Elution Type:** Gradient**Eluent A:** 0.1% TFA/ Water**Eluent B:** 0.085% TFA/ 95% Acetonitrile/Water


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	60	35	65

Flow Rate: 3 µL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Page 1 of 1

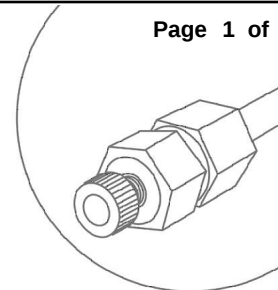

Products used in this Application:

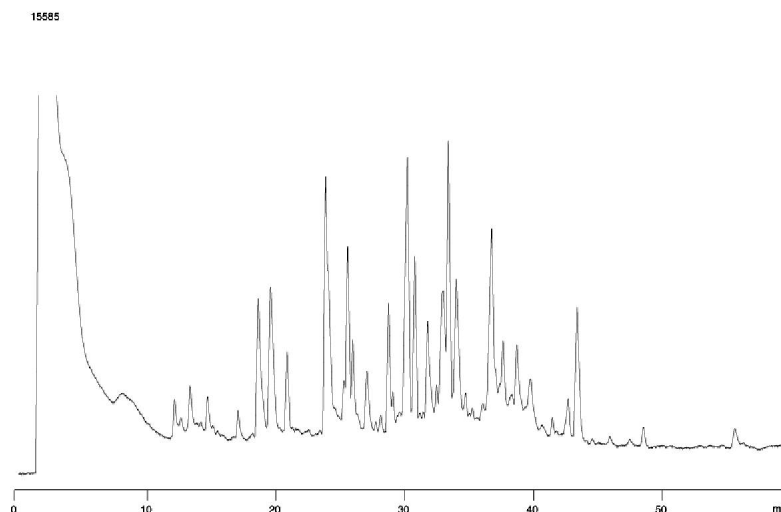ONYX**ANALYTES:**

- 1 Beta Amylase Digest



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.