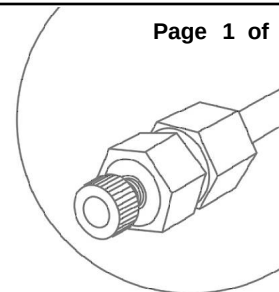
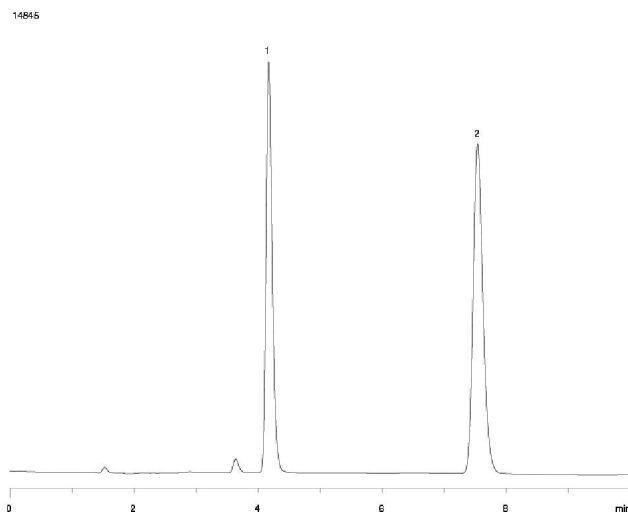


Separation of Fenofibrate and Gemfibrozil using Synergi Fusion-RP

Page 1 of 1

Column: Synergi 4u Fusion-RP 80A**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4424-E0**Elution Type:** Isocratic**Eluent A:** 20mM Potassium Phosphate, pH 2.5**Pct. A:** 30**Eluent B:** Acetonitrile**Pct. B:** 70**Flow Rate:** 1 mL/min**Col. Temp.:** 30 °C**Detection:** UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Products used in this Application:

synergi™

ANALYTES:

- 1 Gemfibrozil
- 2 Fenofibrate



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.