


Abused Inhalants

Column:	Zebron ZB-624	
Dimensions:	30m x 0.32mm x 1.8µm	
Order No.:	7HM-G005-31	
Injection:	Split 20:1 @ 125°C, 100µL headspace	
Carrier Gas:	Constant Flow Helium 3.5 mL/min	
Oven Program:	40°C for 4 min to 120°C at 10°C/min for 4 min.	
Detector:	FID @ 250°C	
Sample:	1. Methanol 2. Ethanol 3. Ether 4. Acetone 5. Isopropanol 6. Methylene chloride 7. Methyl tert butyl ether (MTBE) 8. Hexane 9. Methyl ethyl ketone (MEK) 10. Ethyl acetate 11. Chloroform 12. Benzene 13. Trichloroethylene 14. Methyl Isobutyl Ketone (MIBK) 15. Toluene	

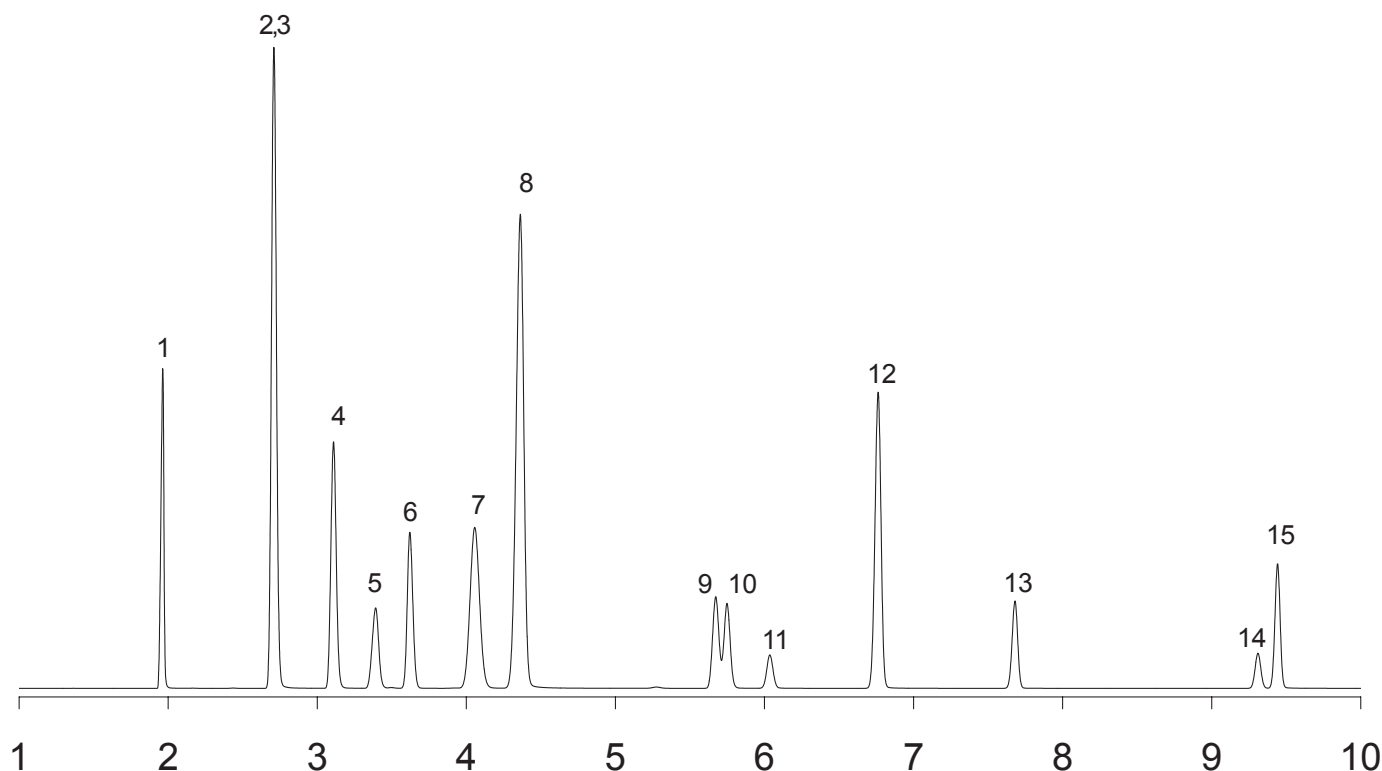

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with tradition™

USA
mail: 411 Madrid Ave.
Torrance, CA
90501-1430
USA

tel: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
273 Sierra Morena,
Suite #104
San Juan,
Puerto Rico 00926

(800) 541-HPLC
(310) 328-7768
info@
phenomenex.com

Canada
411 Madrid Ave.
Torrance, CA
90501-1430
USA

(800) 543-3681
(310) 328-7768
info@
phenomenex.com

France
118/130
avenue Jean Jaurès
75171 Paris Cedex 19
France

01 72 72 92 10
01 72 72 92 11
franceinfo@phenomenex.com

United Kingdom
Queens Avenue,
Hurdfield Ind. Est.,
Macclesfield, Cheshire
SK10 2BN, England

01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
Zeppelinstr. 5
63741 Aschaffenburg
Germany

06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
P O Box 31601
Milford
Auckland
New Zealand

09-4780951
09-4780952
info@
phenomenex.co.nz

Australia
PO Box 4084
Lane Cove, NSW 2066
Australia

02-9428-6444
02-9428-6445
info@
phenomenex.com.au