


Fragrance Analysis

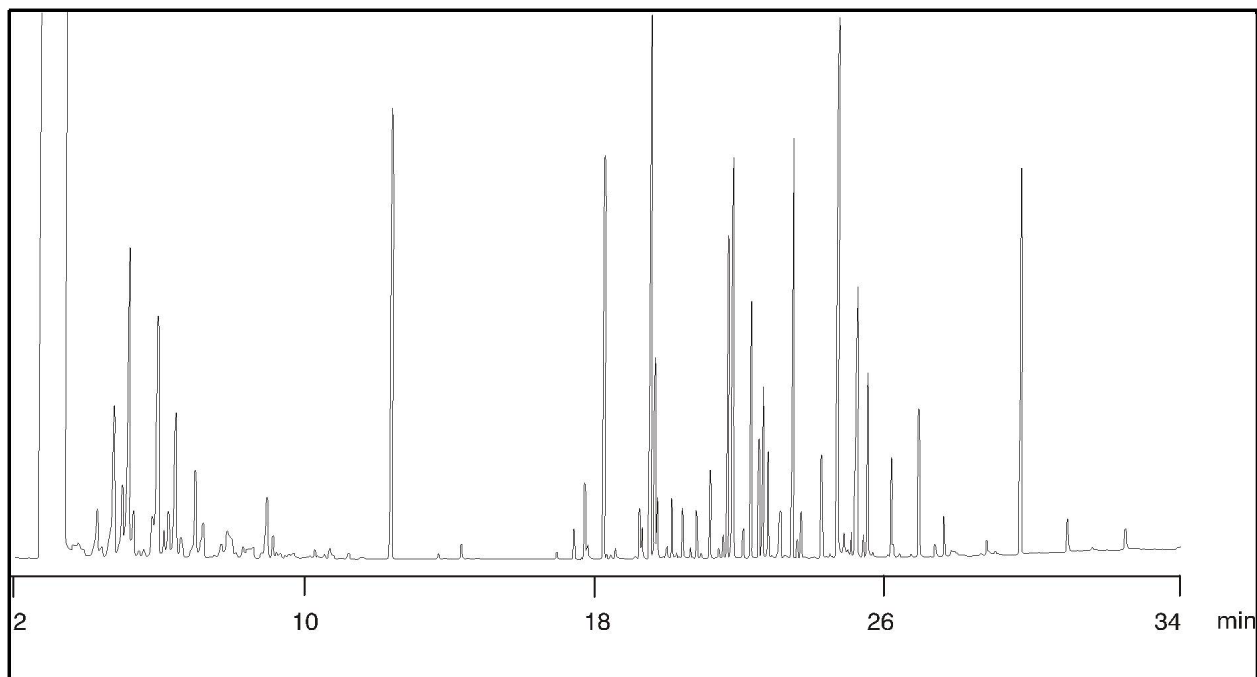

Column: Zebron ZB-WAX
Phase: Polyethylene glycol (PEG)
Dimensions: 30 meters x 0.25 mm x 0.25 μ m
Order No:

Gradient Profile:

Step No.	Step Time	Flow Rate	Temperature
1	0	0	40
2	8	0	40
3	29	0	250
4	31	0	250

Carrier Gas: Helium 30 cm/sec @ 40 °C
Flow Rate: 0.88 mL/min
Injection: Split 50:1 1 μ L @ 250 °C
Detection: Flame Ionization (FID) 300 °C

SAMPLE ANALYTE

1 Unidentified

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1308/ZEB/013 / GC_17a.wmf