


APPLICATION NOTE - HPLC
IgG Dog Native Acetonitrile



14918

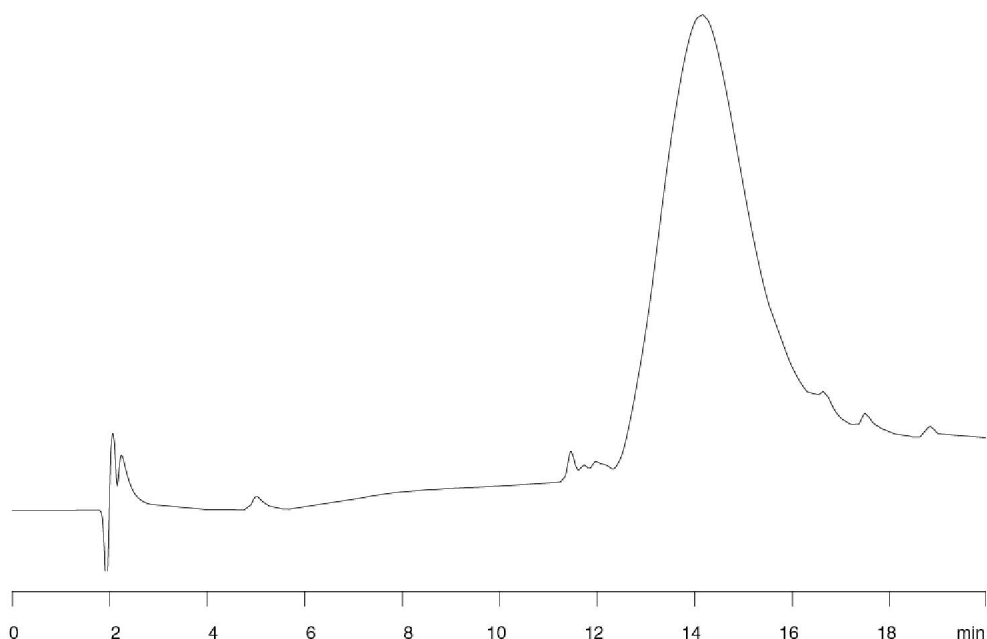

Column: Jupiter 5 μ m C4
Dimensions: 150 x 2 mm ID
Order No: 00F-4167-B0
Elution Type: Gradient
Eluent A: 0.1% TFA in 95:5 Water/ACN
Eluent B: 0.085% TFA in 95:5 ACN/Water

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	80	20
2	20	5	95

Flow Rate: 0.25 mL/min
Temperature: 22 °C
Detection: UV-Vis Abs. - Diode Array @ 220 nm
SPE:

SAMPLE ANALYTES:

1 IgG_Dog Native

For more information call your Phenomenex Representative at (310) 212-0555

©Copyright 2000-2004 Phenomenex Inc.

Ref. No.: 14918

Page 1 of 1

