


Odd and Even Carbon FAMES

Column:	Zebron ZB-5	
Dimensions:	30m x 0.25mm x 0.25µm	
Order No.:	7HG-G002-11	
Injection:	Split 20:1 @ 280°C, 0.2µL	
Carrier Gas:	Constant Flow Helium 1.1 mL/min	
Oven Program:	100°C for 2min to 280°C at 8°C/min for 5 min	
Detector:	MS @ 280°C; 25-400amu	
Sample:	1. Methyl Caproate	11. Methyl Heptadecanoate
	2. Methyl Caprylate	12. Methyl Stearate
	3. Methyl Nonoate	13. Methyl Nonadecanoate
	4. Methyl Caprate	14. Methyl Arachidate
	5. Methyl Undecanoate	15. Methyl Heneicosanoate
	6. Methyl Laurate	16. Methyl Behenate
	7. Methyl Tridecanoate	17. Methyl Tricosanoate
	8. Methyl Myristate	18. Methyl Lignocerate
	9. Methyl Pentadecanoate	
	10. Methyl Palmitate	

Notes: All analytes were 2.8mg/mL except 18 at 0.7mg/mL.

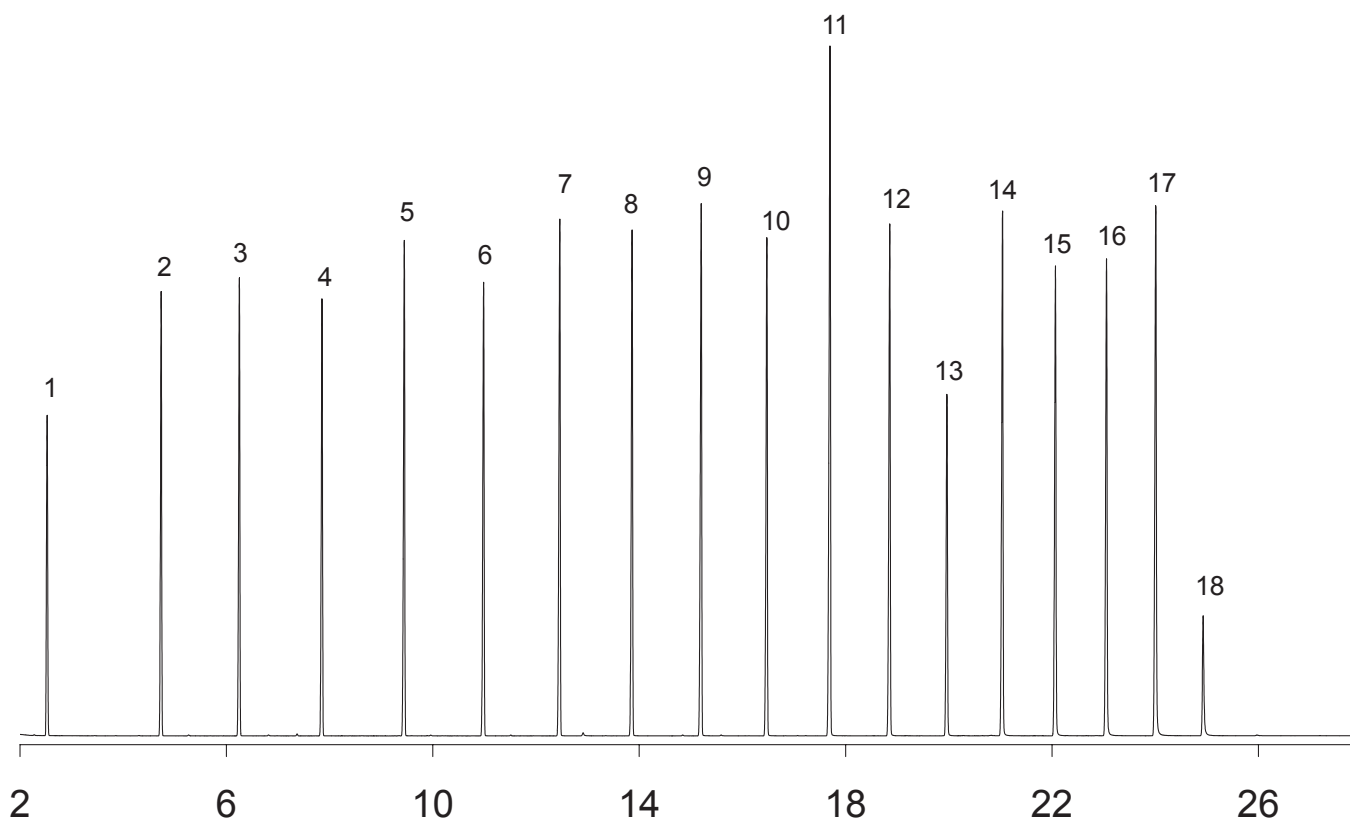

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with traditionSM

USA mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	Puerto Rico 271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926 (800) 541-HPLC (310) 328-7768 info@phenomenex.com	Canada 411 Madrid Ave. Torrance, CA 90501-1430 USA (800) 543-3681 (310) 328-7768 info@phenomenex.com	United Kingdom Queens Avenue, Hurdfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England 01625-501367 01625-501796 ukinfo@phenomenex.com	Germany Zeppelinstr. 5 63741 Aschaffenburg Deutschland 06021-58830-0 06021-58830-11 anfrage@phenomenex.com	New Zealand P.O. Box 31-601 Milford Auckland New Zealand 09-4780951 09-4780952 info@phenomenex.co.nz	Australia P.O. Box 264 Pennant Hills NSW 1715 Australia 1800-553-929 1800-553-923 info@phenomenex.com.au
---	---	--	---	---	--	--