

IgG Dog Rd_2

Column: Jupiter 5u C4 300A

Dimensions: 150 x 2 mm ID

Order No: 00F-4167-B0

Elution Type: Gradient

Eluent A: 0.1% TFA in 95:5 Water/ACN


Eluent B: 0.085% TFA in 75:20:5 ACN/IPA/Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	20	10	90


Flow Rate: 0.25 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Diode Array @ 220 nm (ambient)

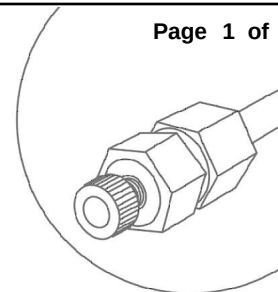

Products used in this Application:

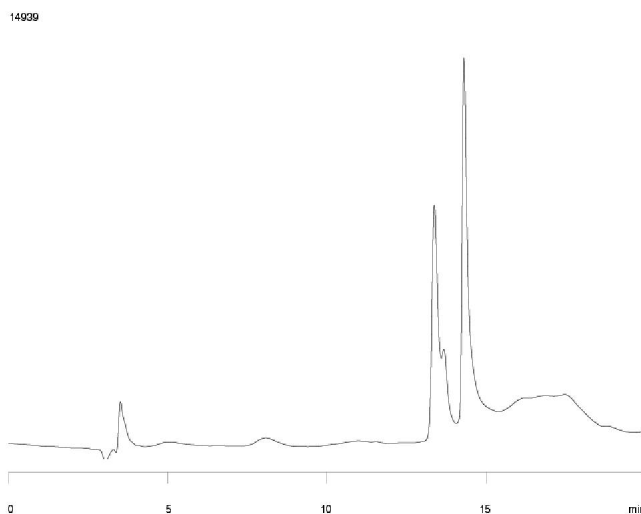

ANALYTES:

- 1 IgG Dog Reduced



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.