


Amphetamines (PFPA Derivatization)

Column: Zebron ZB-1ms
Dimensions: 20meter x 0.18mm x 0.18µm
Order No.: 7FD-G011-08
Injection: Split 10:1 @ 290°C, 1µL
Carrier Gas: Constant Flow Helium 0.8mL/min
Oven Program: 80°C to 300°C at 20°C/min
Detector: MSD; 35-500 amu
Sample: 1. Amphetamine (PFP)
 2. Methamphetamine (PFP)

Gold Seal: Phenomenex single groove (AG0-7518)
Liner: Split/Splitless Single Taper (AG0-7516)
Septum: Phenomenex PhenoRed-400 (AG0-4690)

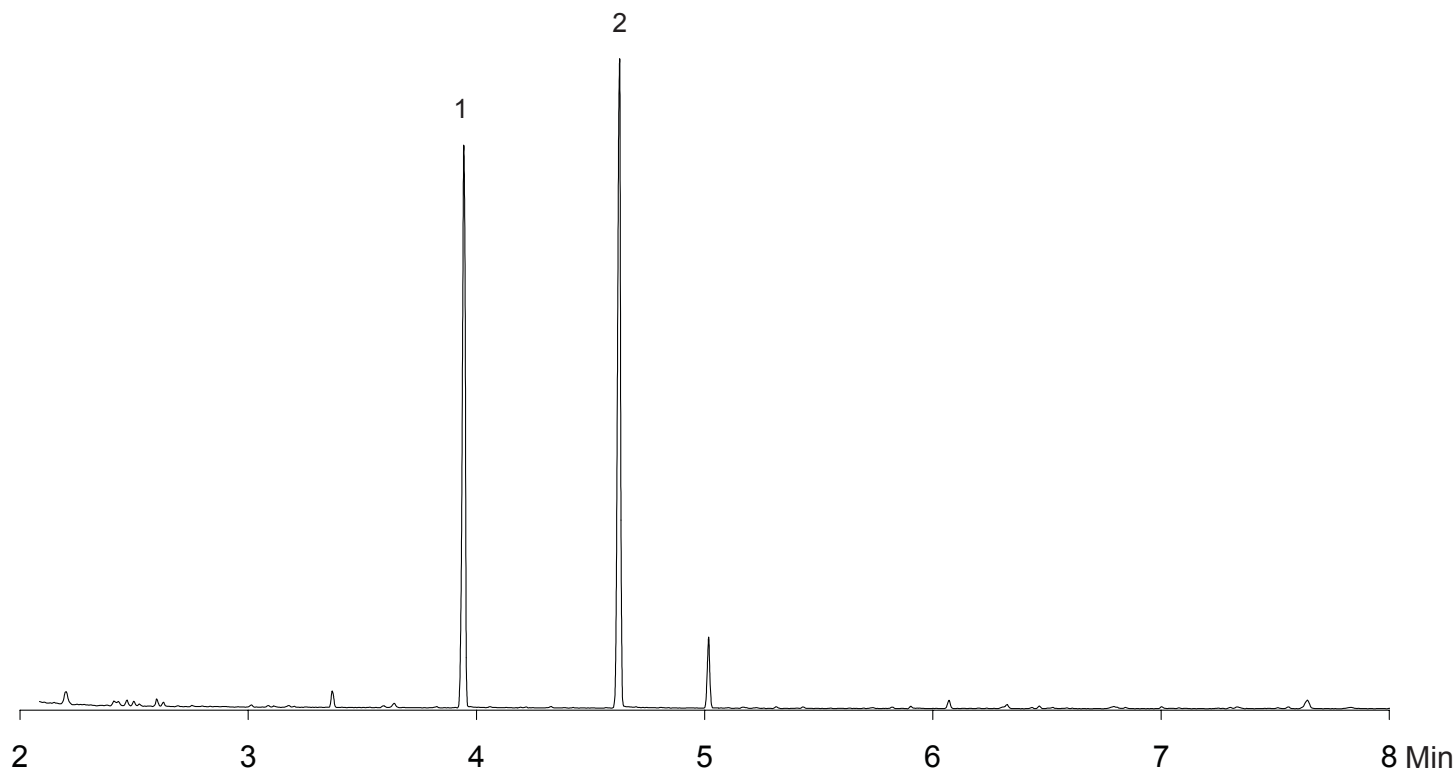

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.



...breaking with tradition™

USA	Puerto Rico	Canada	United Kingdom	Germany	New Zealand	Australia
mail: 411 Madrid Avenue Torrance, CA 90501-1430 USA	273 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926	411 Madrid Avenue Torrance, CA 90501-1430 USA	Queens Avenue, Hurdfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England	Zeppelinstr. 5 63741 Aschaffenburg Deutschland	P.O. Box 31-601 Milford Auckland New Zealand	P.O. Box 4084 Lane Cove, NSW 1595 Australia
tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	(800) 541-HPLC (310) 328-7768 info@phenomenex.com	(800) 543-3681 (310) 328-7768 info@phenomenex.com	01625-501367 01625-501796 ukinfo@phenomenex.com	06021-58830-0 06021-58830-11 anfrage@phenomenex.com	09-4780951 09-4780952 info@phenomenex.co.nz	02-9428-6444 02-9428-6445 info@phenomenex.com.au