


APPLICATION NOTE - HPLC
Cephalosporin Antibiotics



14997

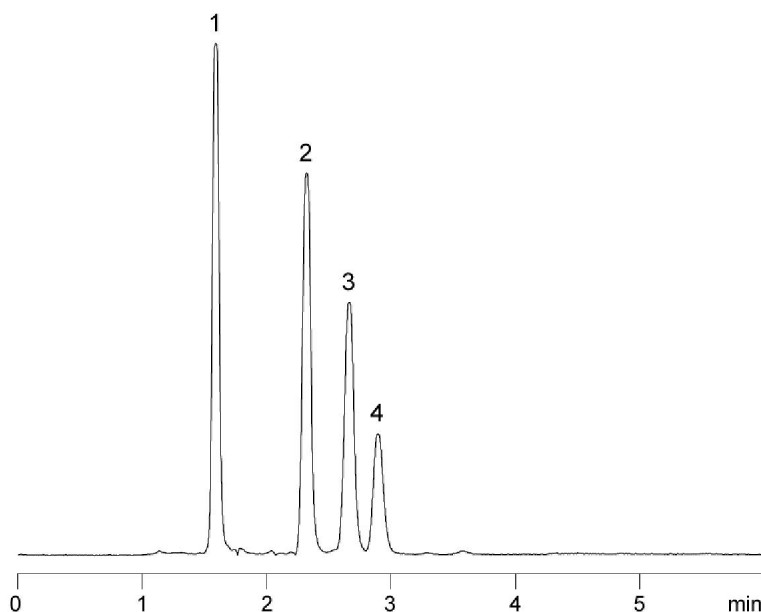

Column: Gemini 5 μ m C18 / ODS C18

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic

Eluent A: 0.5% Acetic acid

Pct. A: 85

Eluent B: Acetonitrile

Pct. B: 15

Flow Rate: 1 mL/min

Temperature: 25 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTES:

- 1 Cefadroxil
- 2 Cephalexin
- 3 Cephadrine
- 4 Cefaclor

SPE:

For more information call your Phenomenex Representative at

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 14997/BRO

Page 1 of 1

