


APPLICATION NOTE - HPLC
Local Anesthetics



14999

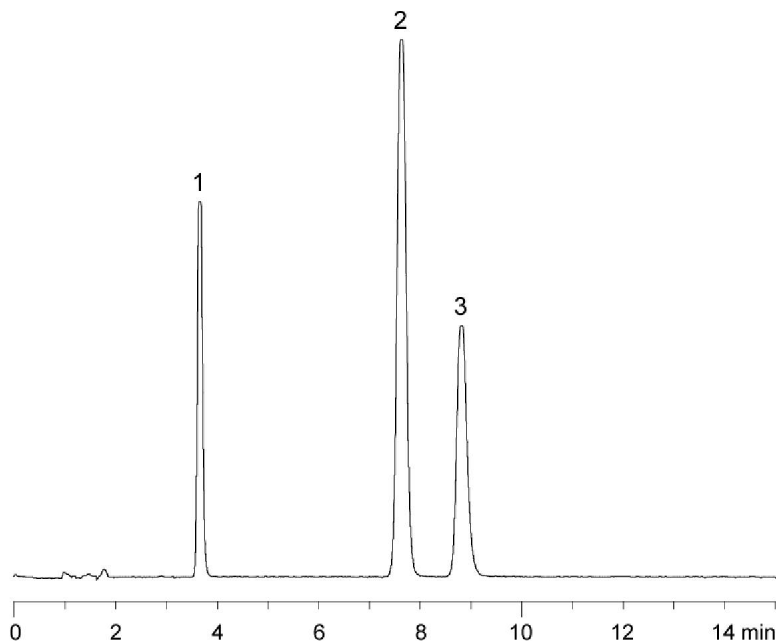

Column: Gemini 5 μ m C18 / ODS C18

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic

Eluent A: 10mM NH_4HCO_3 , pH 10.5

Pct. A: 50

Eluent B: Acetonitrile

Pct. B: 50

Flow Rate: 1 mL/min

Temperature: 25 $^{\circ}\text{C}$

Detection: UV-Vis Abs. - Variable Wave. @ 210 nm

SAMPLE ANALYTES:

- 1 Procaine
- 2 Lidocaine
- 3 Tetracaine

SPE:

For more information call your Phenomenex Representative at

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 14999/BRO

Page 1 of 1

