

White wine analysis using Onyx C18

Page 1 of 1

Column: Onyx Monolithic C18

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7643

Elution Type: Gradient

Eluent A: 0.1% H₃PO₄ in Water

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	10	97	3
	2	22	82	18
	4	30	50	50
	5	35	50	50

Flow Rate: 2 mL/min


Col. Temp.: 30 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Products used in this Application:

ONYX

15591

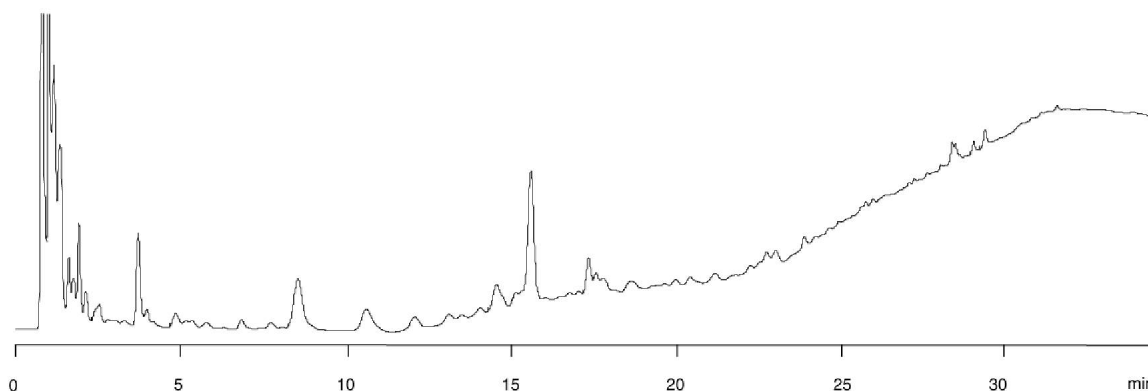

ANALYTES:

- 1 White wine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.