


Barbiturates and Cocaine

Column: Zebron ZB-5
Dimensions: 30meter x 0.25mm x 0.25µm
Order No.: 7HG-G002-11
Injection: Splitless @ 250°C, 2µL
Carrier Gas: Constant Flow Helium 3.0 mL/min
Oven Program: 75°C to 320°C at 25°C/min.
Detector: MSD; scan 50-500
Sample:

1. Pentobarbital
2. Meprobamate
3. Hexobarbital
4. Cocaine
5. Diazepam
6. 6-MAM

Notes: Analytes are 1µg/100µL. Analytes were extracted from serum using Strata™ Screen-C GF and derivatized to form -TMS derivatives.

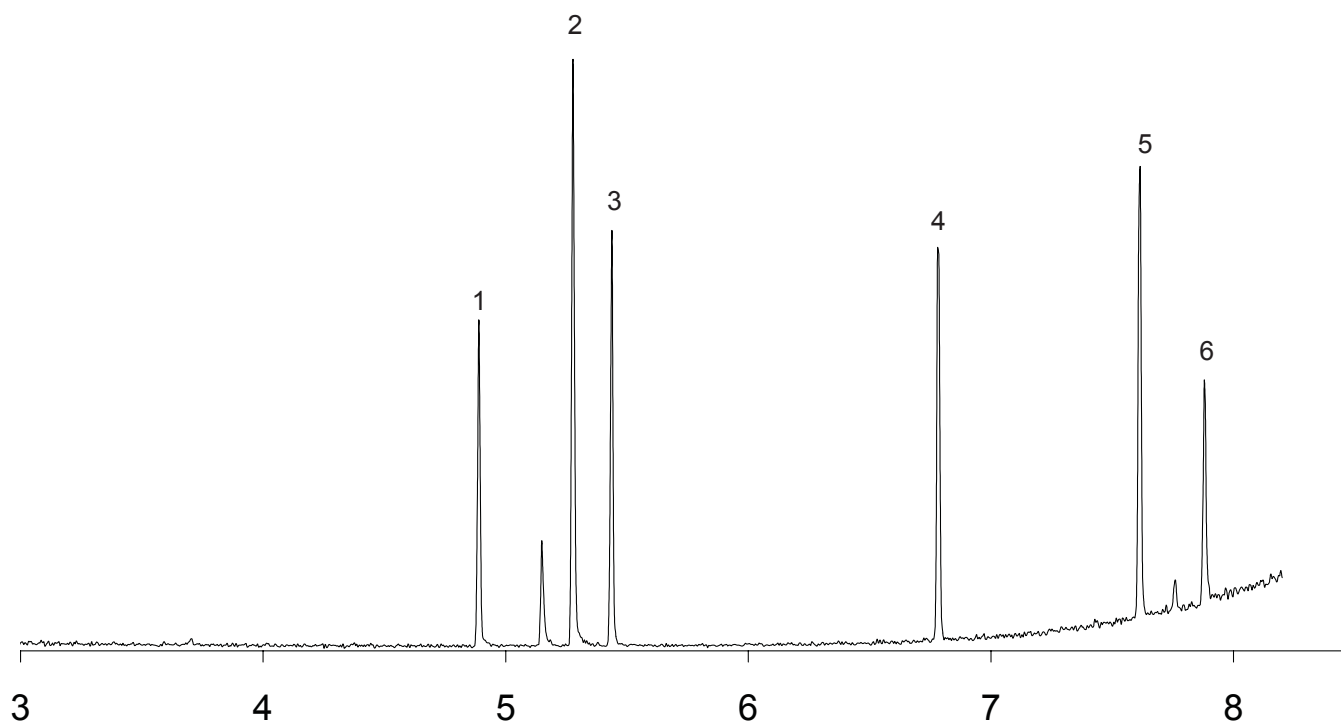

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.


phenomenex®
 ...breaking with traditionSM

USA	Puerto Rico	Canada	United Kingdom	Germany	New Zealand	Australia
mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926 (800) 541-HPLC (310) 328-7768 info@phenomenex.com	411 Madrid Ave. Torrance, CA 90501-1430 USA (800) 543-3681 (310) 328-7768 info@phenomenex.com	Queens Avenue, Hurdfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England 01625-501367 01625-501796 ukinfo@phenomenex.com	Zeppelinstr. 5 63741 Aschaffenburg Deutschland 06021-58830-0 06021-58830-11 anfrage@phenomenex.com	P.O. Box 31-601 Milford Auckland New Zealand 09-4780951 09-4780952 info@phenomenex.co.nz	P.O. Box 264 Pennant Hills NSW 1715 Australia 1800-553-929 1800-553-923 info@phenomenex.com.au