

'Soy Complex' Extract (Finished Product)

Page 1 of 1

Column: Onyx Monolithic C18

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7643

Elution Type: Gradient

Eluent A: 0.1% TFA in Water

Eluent B: 0.08% TFA in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	15	30	70

Flow Rate: 1 mL/min

Col. Temp.: 30 °C


Detection: UV-Vis Abs. - Variable Wave. @ 280 nm (ambient)

Products used in this Application:

ONYX[™]

Analyst Note: 2 columns coupled end to end

15593

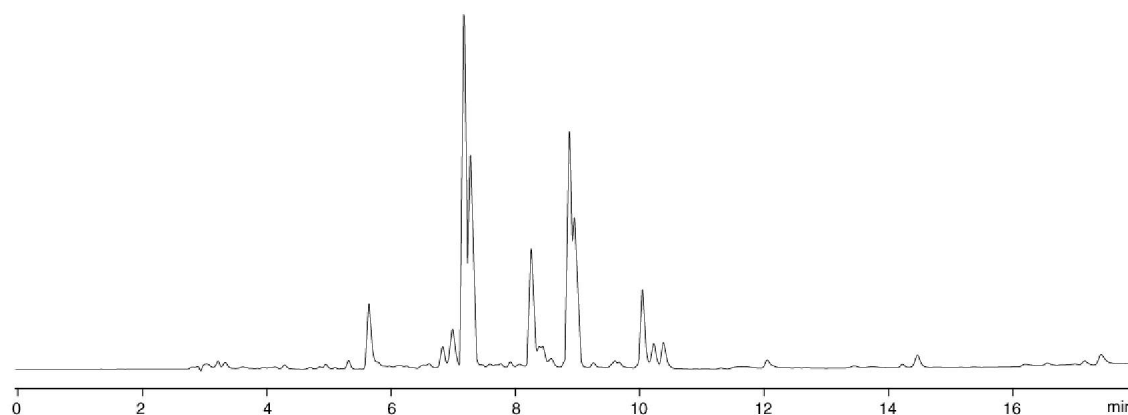

ANALYTES:

1 Soy Complex



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.