


Brominated Diphenyl Ethers (BDE's)

Column:	Zebron ZB-5		
Dimensions:	30m x 0.25mm x 0.25µm		
Order No.:	7HG-G002-11		
Injection:	Splitless pulsed 30psi @ 275°C, 1µL		
Carrier Gas:	Programmed Flow Helium 1mL/min for 45.74 min to 2.0mL/min at 70mL/min		
Oven Program:	110°C for 1 min to 180°C at 12°C/min for 1 min to 240°C at 4°C/min for 5min to 280°C at 4°C/min for 10 min.		
Detector:	MS; 50-800 amu		
Sample:	1. BDE 3	11. BDE 66	21. BDE 153
	2. BDE 7	12. Mirex	22. HBCD
	3. BDE 15	13. PCB 198	23. BDE 138
	4. PCB 112	14. BDE 77	24. BDE 166
	5. BDE 17	15. BDE 100	25. BDE 128
	6. BDE 28	16. BDE 119	26. BDE 183
	7. BDE 50	17. BDE 99	27. BDE 181
	8. BDE 49	18. BDE 85	
	9. BDE 71	19. BDE 126	
	10. BDE 47	20. BDE 154	

Notes: BDE Wellington Mix at 1µg/mL

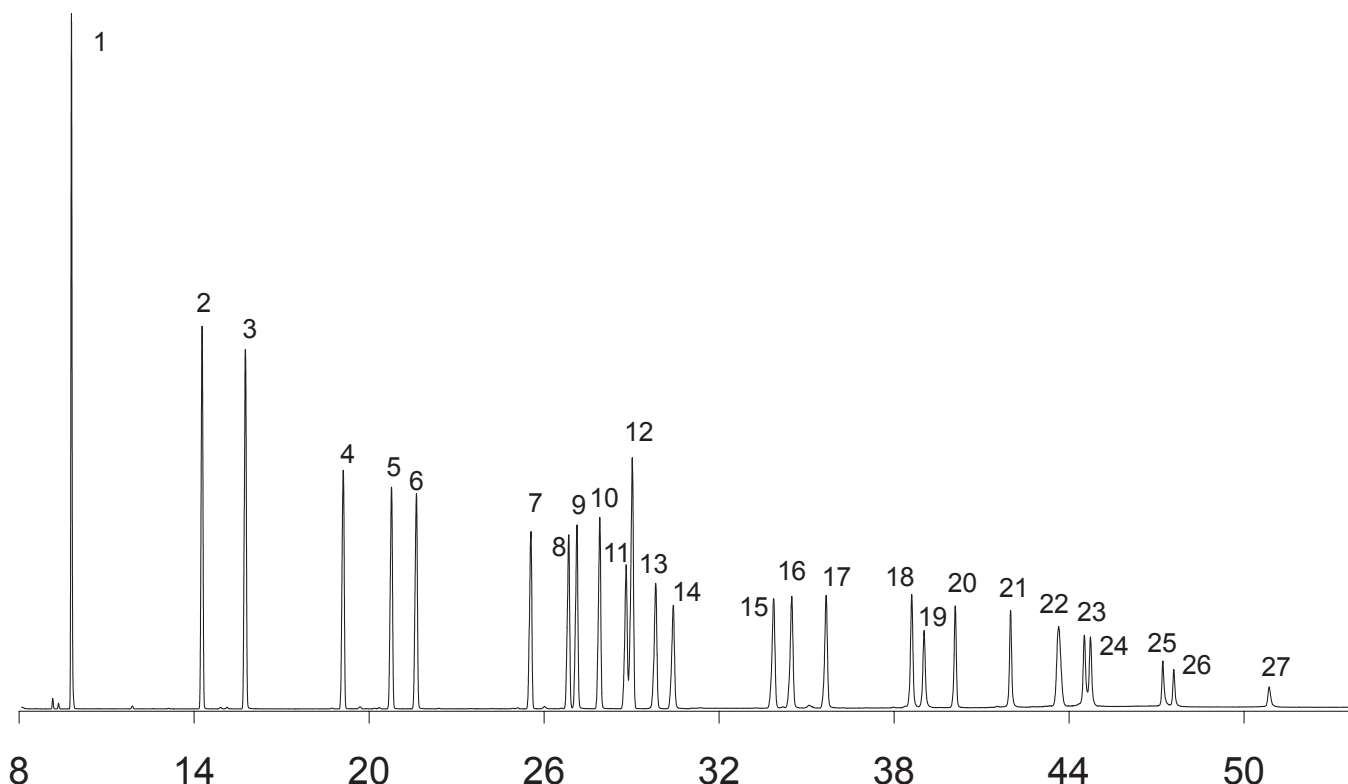

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with traditionSM

USA mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	Puerto Rico 271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926 (800) 541-HPLC (310) 328-7768 info@phenomenex.com	Canada 411 Madrid Ave. Torrance, CA 90501-1430 USA (800) 543-3681 (310) 328-7768 info@phenomenex.com	United Kingdom Queens Avenue, Hurdsfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England 01625-501367 01625-501796 ukinfo@phenomenex.com	Germany Zeppelinstr. 5 63741 Aschaffenburg Deutschland 06021-58830-0 06021-58830-11 anfrage@phenomenex.com	New Zealand P.O. Box 31-601 Milford Auckland New Zealand 09-4780951 09-4780952 info@phenomenex.co.nz	Australia P.O. Box 264 Pennant Hills NSW 1715 Australia 1800-553-929 1800-553-923 info@phenomenex.com.au
--	--	---	---	---	---	--