


Cough and Cold

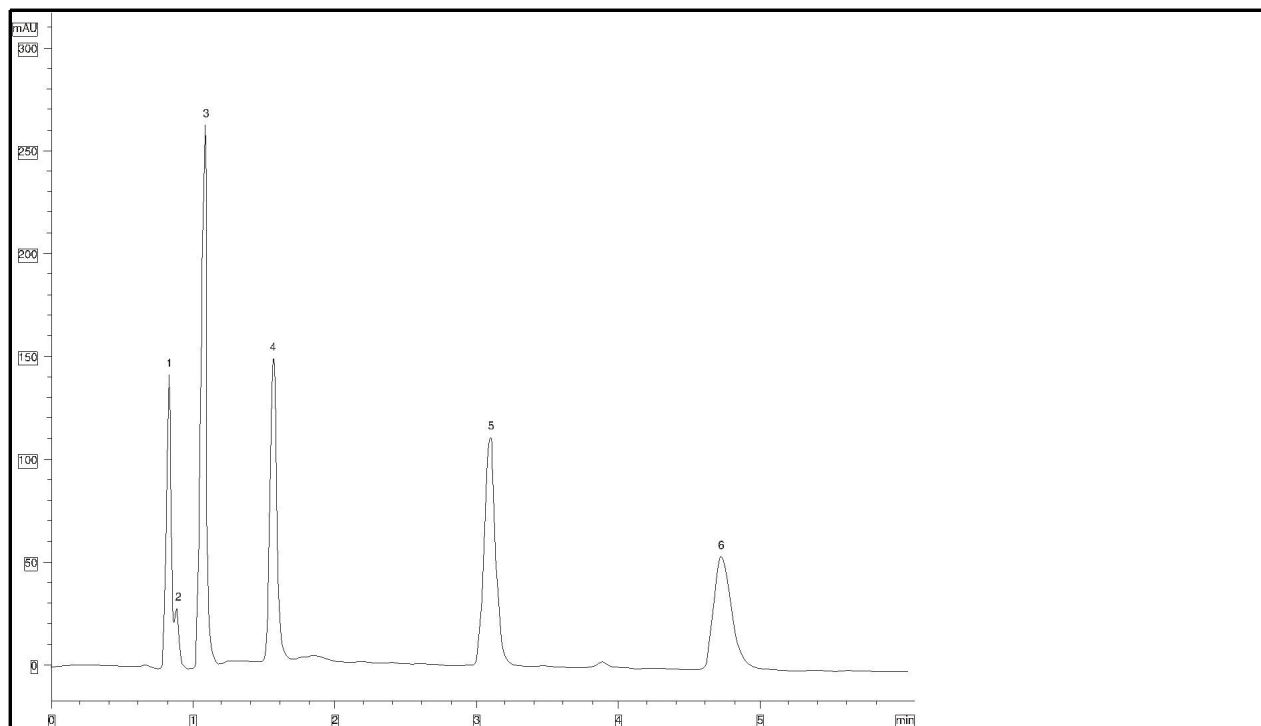

Column: Luna 5 μ m C6H5 / Phenyl Hexyl

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 20 mM Potassium phosphate pH 10.0

Pct. A: 35

Eluent B: Acetonitrile

Pct. B: 35

Eluent C: Methanol

Pct. C: 30

Flow Rate: 1.5 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 214 nm

SAMPLE ANALYTE

- 1 HBr
- 2 Maleic acid
- 3 Acetaminophen
- 4 Pseudoephedrine
- 5 Chlorpheniramine
- 6 Dextromethorphan

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 2574/ALG/184