

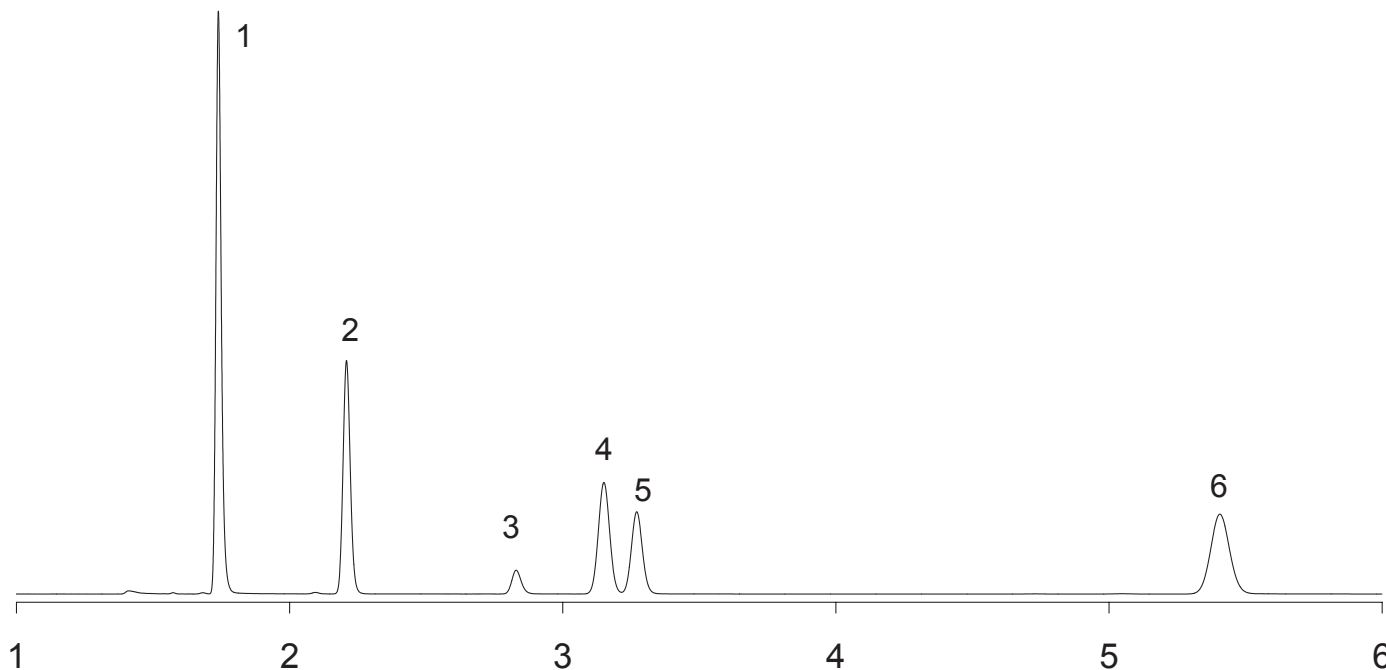
Blood Alcohols

App ID

Column: Zebtron ZB-WAX
Dimensions: 30m x 0.53mm x 1.0µm
Order No.: 7HK-G007-22
Injection: Splitless Headspace @ 225°C, 1mL
Carrier Gas: Helium 4.7mL/min
Oven Program: 55°C for 6 min.
Detector: FID @ 325°C
Sample:

1. Acetaldehyde
2. Acetone
3. Methanol
4. Isopropanol
5. Ethanol
6. n-Propanol

Notes: 1 mL of headspace was taken from sample heated for 10 minutes at 70°C.


www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.


phenomenex
...breaking with traditionSM

USA	Puerto Rico	Canada	United Kingdom	Germany	New Zealand	Australia
mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA	271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926	411 Madrid Ave. Torrance, CA 90501-1430 USA	Queens Avenue, Hurdsfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England	Zeppelinstr. 5 63741 Aschaffenburg Deutschland	P.O. Box 31-601 Milford Auckland New Zealand	P.O. Box 264 Pennant Hills NSW 1715 Australia
tel.: (310) 212-0555	(800) 541-HPLC	(800) 543-3681	01625-501367	06021-58830-0	09-4780951	1800-553-929
fax: (310) 328-7768	(310) 328-7768	(310) 328-7768	01625-501796	06021-58830-11	09-4780952	1800-553-923
email: info@phenomenex.com	info@phenomenex.com	info@phenomenex.com	ukinfo@phenomenex.com	anfrage@phenomenex.com	info@phenomenex.co.nz	info@phenomenex.com.au