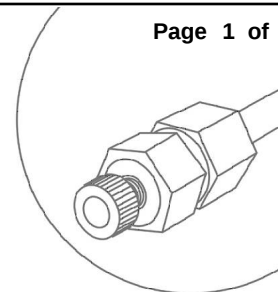
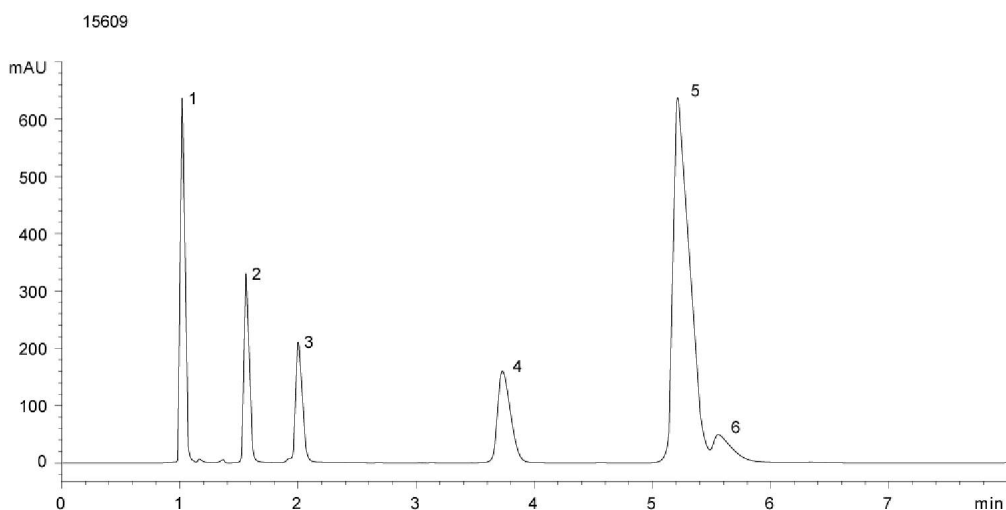


Beta-Blockers (polar basics) at low pH by Gemini 5um C18

Page 1 of 1

Column: Gemini 5u C18 110A**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4435-E0**Elution Type:** Isocratic**Eluent A:** 20mM Phosphate buffer pH2.5**Pct. A:** 75**Eluent B:** Acetonitrile**Pct. B:** 25**Flow Rate:** 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 230 nm (ambient)**Products used in this Application:****Analyst Note:** Peak 4 -- Labetolol
Efficiency = 7127 Peak Symmetry = 0.68**ANALYTES:**

- 1 Atenolol
- 2 Pindolol
- 3 Metoprolol
- 4 Labetolol
- 5 Propranolol
- 6 Alprenolol

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country,
contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.