

Cough And Cold USP Method

Column: Luna 5u Phenyl-Hexyl 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4257-E0

Elution Type: Isocratic

Eluent A: Methanol

Eluent B: Water w/0.1% H₃PO₄ & 0.1% heptane sulfonate

Eluent C: Acetonitrile w/0.1% H₃PO₄

Flow Rate: 2.05 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 214 nm (ambient)

Pct. A: 35

Pct. B: 55


Pct. C: 10

Page 1 of 1

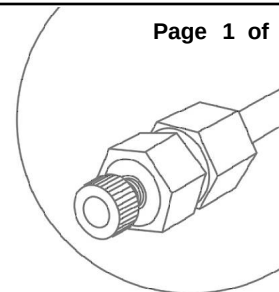

Products used in this Application:

LUNA™

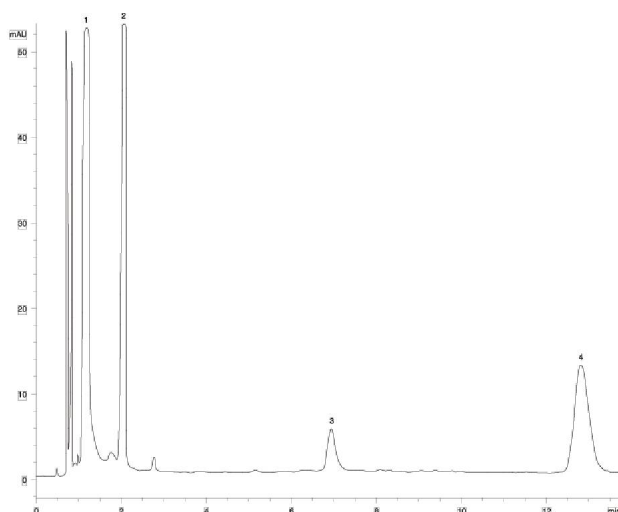

ANALYTES:

- 1 Acetaminophen
- 2 Pseudoephedrine
- 3 Chlorpheniramine
- 4 Dextromethorphan

